

ENT

Endolaryngeal Microsurgery



DISTRIBUTED BY:



SURGICAL ASIA Sdn Bhd (1272656-H)

Suppliers of Fine Quality Instruments

No. 8, Lebuhr Perusahaan Klebang 9A,
Kawasan Perusahaan IGB, 31200 Chemor, Perak, Malaysia

Phone: +605 210 8080 • www.inkasurgical.asia



Exclusive Distributor

OBEX Medical Ltd

0800 656 239

www.inkasurgical.co.nz



we thank you for your support

Disclaimer: Product availability is subject to regulatory approval in individual markets.
Not all products are available in all markets. Contact us to find out more.

BRC-00073

contents

01	Legend	25	Needleholder - Vocal Cord Spreader
03	Chest Supports	26	Foreign Body Forceps
04	Operating Laryngoscopes	28	Esophagoscopy Forceps
07	Modified Laryngoscopes	30	Hand Instruments
08	Laryngoscopes	31	Suction Tubes
11	Cylindrical Shafts	32	Mirrors
14	Conical Shafts	34	Laryngeal Applicators
15	Reverse Action Handles	35	Scissors - Forceps - Retractors
16	Paediatric Patterns	36	SUTTER Monopolar Forceps
19	Miniature Range	37	SUTTER Calvin® Bipolar Forceps
21	Cleaning Channel Range	39	SUTTER Curis® RF - Larynx - Technology
23	Cup Forceps - Scissors		

legend

Triple Connector



ACMI




WOLF




STORZ
OLYMPUS

British standard **Triple Connector** with adapters to fit light cables for fiber optic light carriers


Distal End

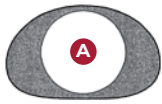


Size 1 -
Infant Size

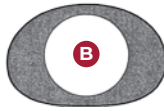


Size 2 -
Child Size

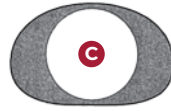


Size 3 -
Adult, Small

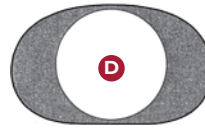


Size 4 -
Adult, Medium

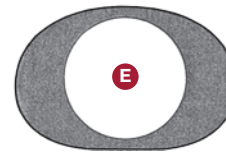

Size 5 -
Adult, Large

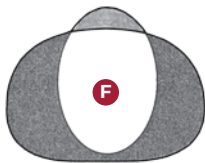


Size 6 -
Anterior
Commissure

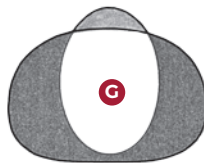


Size 7 -
for Long Neck

Distal ends have 10° upwards curve to allow better access to anterior commissure.

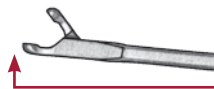

Shaft

Shaft refers to the long cylindrical part of an instrument. Example on the right.

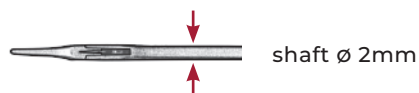

legend

Reverse Action Handles

Reverse Action Handles are thumb activated and reduces distal tip movement.



FG Sizes

Figure (FG) sizes and their corresponding diameter sizes.

FG	Ø mm
6	Ø 2
7	Ø 2.3
8	Ø 2.7
9	Ø 3
10	Ø 3.3
12	Ø 4
15	Ø 5

Laryngeal & Rhinoscopic Mirrors



Size 0 -
Ø 10mm




Size 1 -
Ø 12mm




Size 2 -
Ø 14mm



Size 3 -
Ø 16mm

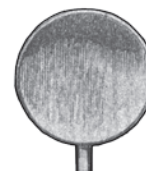


Size 4 -
Ø 18mm

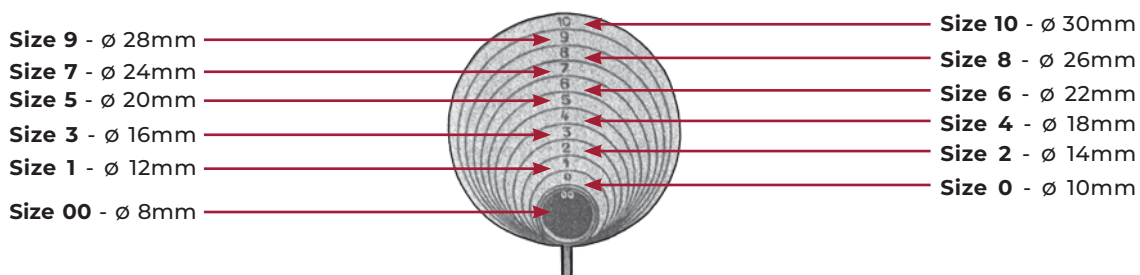

chest supports

Instruments @ 50% · Inserts to Scale

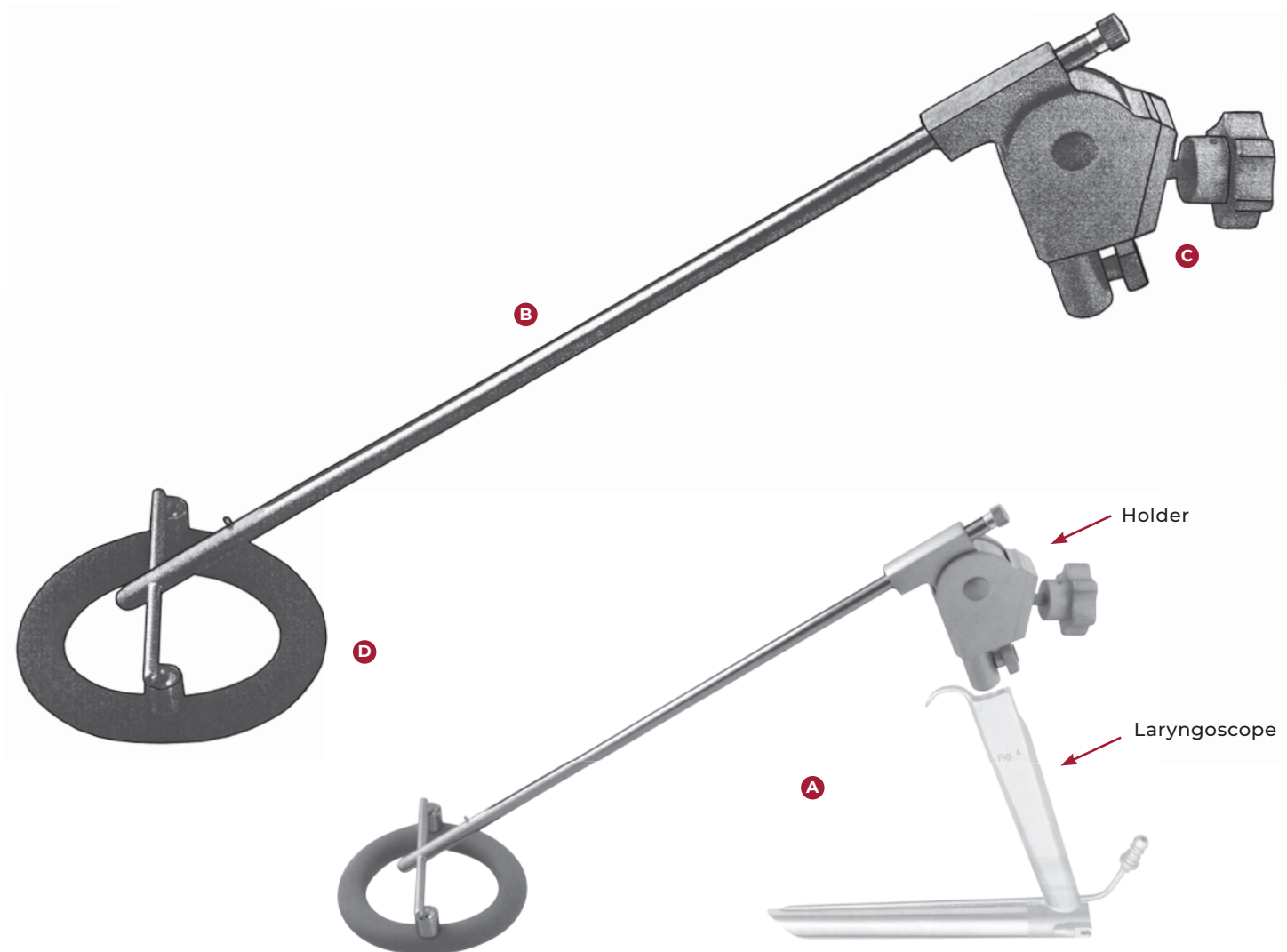

A 35400.34
RIECKER-KLEINSASSER
 Laryngoscope
 with rod and holder, 34cm rod length

C 35400.50
 Laryngoscope Holder Only
 for chest support

B 35403.24
 Rod only with Rubber Ring
 for children, used with **35400.34**
 24cm length


D 35400.99
 Spare Rubber Ring
 latex free for chest support

operating laryngoscopes


Instruments @ 50% · Inserts to Scale · fiber optics not included



A 35390.01
KLEINSASSER style
Size 1, infant size
I.D. 11mm - 135mm

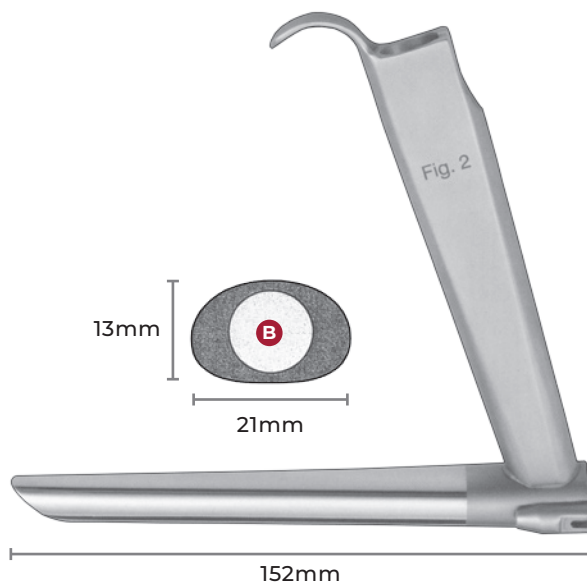


B 35390.02
KLEINSASSER style
Size 2, child size
I.D. 11mm - 152mm

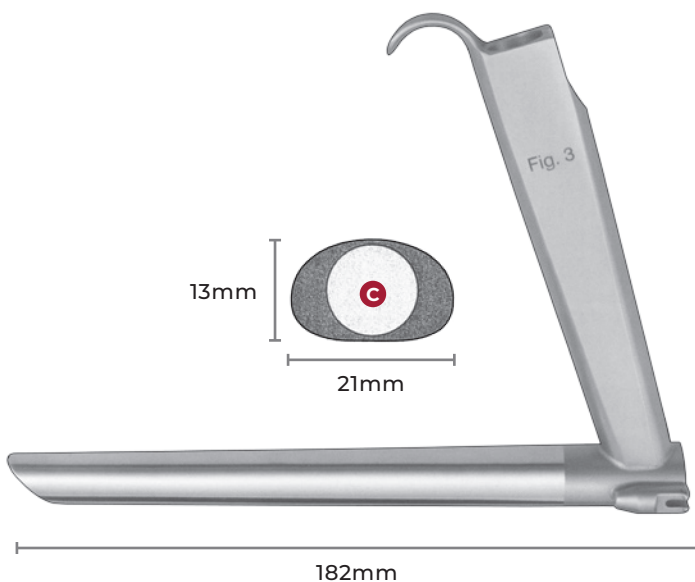


C 35390.03
KLEINSASSER style
Size 3, adult size, small
I.D. 12mm - 182mm

operating laryngoscopes

Instruments @ 50% · Inserts to Scale · fiber optics not included

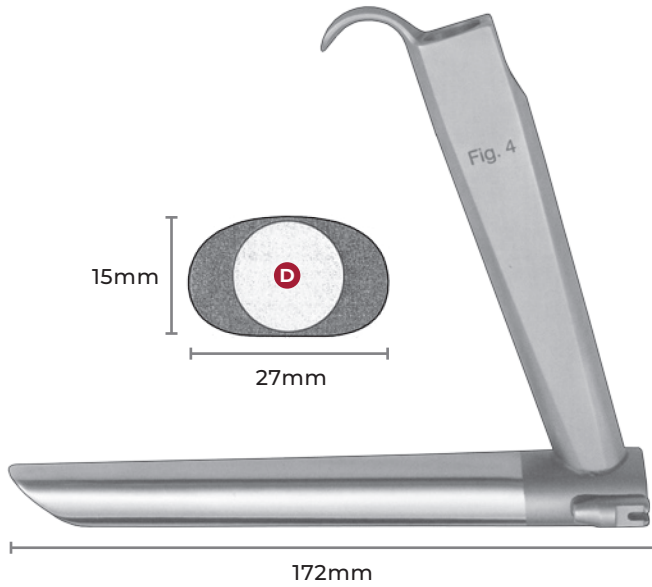


- D 35390.04**
KLEINSASSER style
Size 4, adult size, medium
I.D. 14mm - 172mm

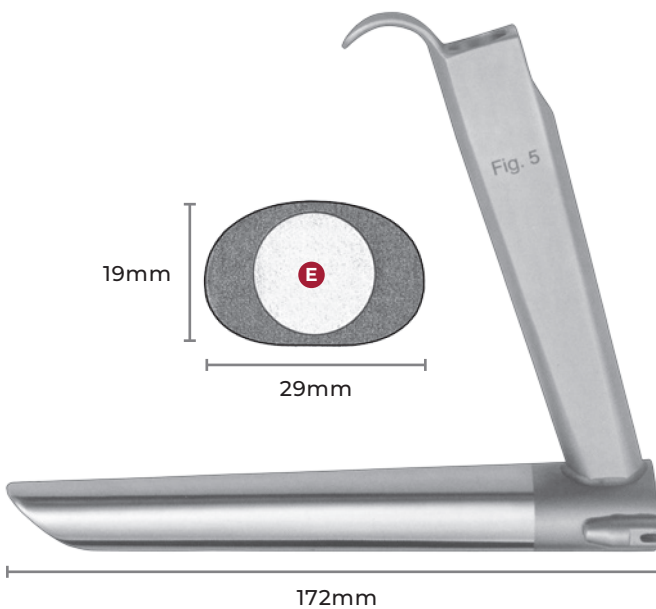


- E 35390.05**
KLEINSASSER style
Size 5, adult size, large
I.D. 16mm - 172mm

operating laryngoscopes

Instruments @ 50% · Inserts to Scale · fiber optics not included

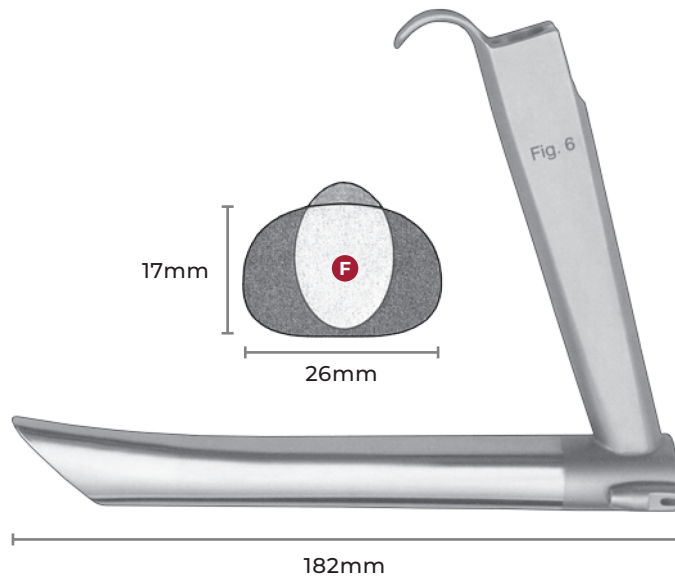


- F 35390.06**
KLEINSASSER style
Size 6, anterior commissure
I.D. 13/18mm - 182mm

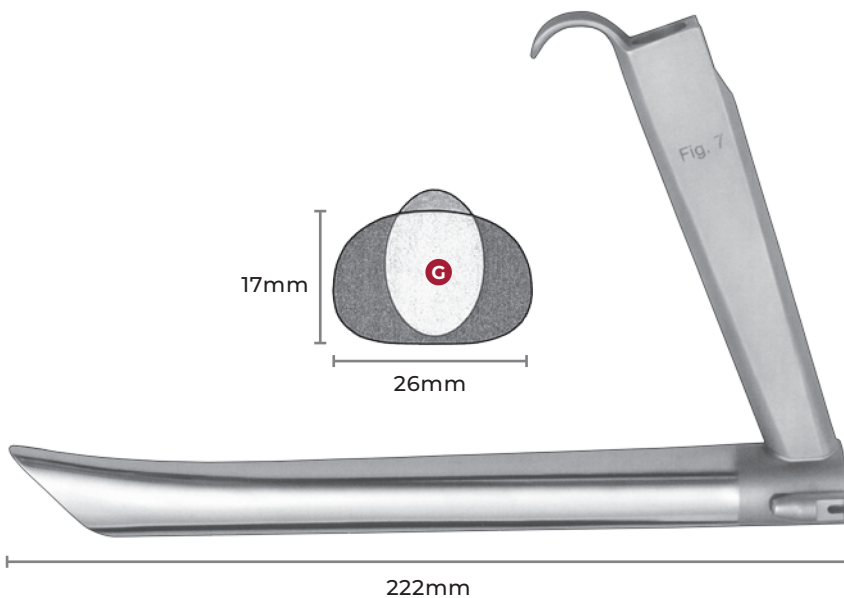


- G 35390.07**
KLEINSASSER style
Size 7, for long necks
I.D. 13/18mm - 222mm

modified laryngoscopes

Instruments @ 50% · Inserts to Scale

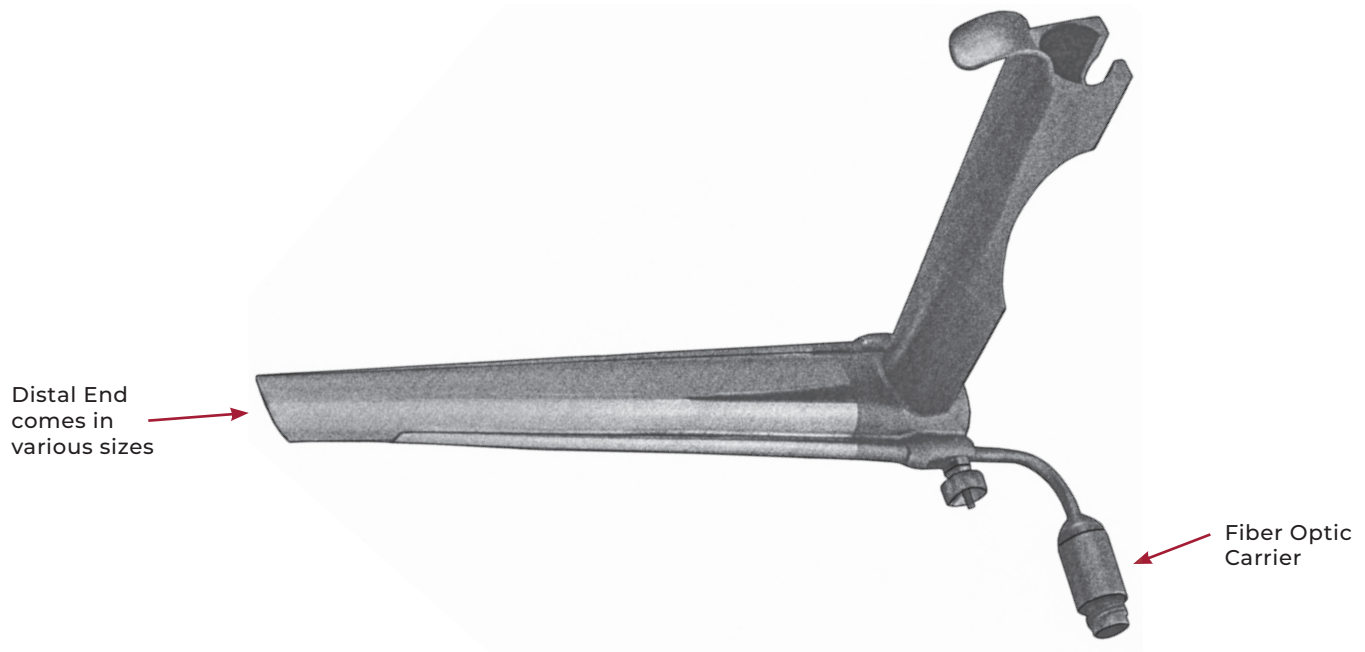

Modified for CO₂ Laser Surgery

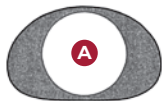

35392.01
KLEINSASSER style
Size 1, infant size
I.D. 11mm - 135mm




35392.02
KLEINSASSER style
Size 2, child size
I.D. 11mm - 152mm




35392.03
KLEINSASSER style
Size 3, adult size, small
I.D. 12mm - 182mm



35392.04
KLEINSASSER style
Size 4, adult size, medium
I.D. 14mm - 172mm

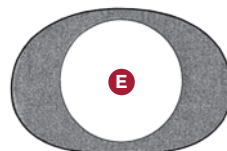

35392.05
KLEINSASSER style
Size 5, adult size, large
I.D. 16mm - 172mm

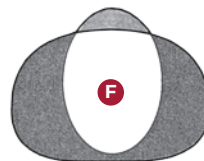


35392.06
KLEINSASSER style
Size 6, anterior commissure
I.D. 13/18mm - 182mm




35392.07
KLEINSASSER style
Size 7, for long necks
I.D. 13/18mm - 222mm

laryngoscopes

Instruments @ 50% · Inserts to Scale

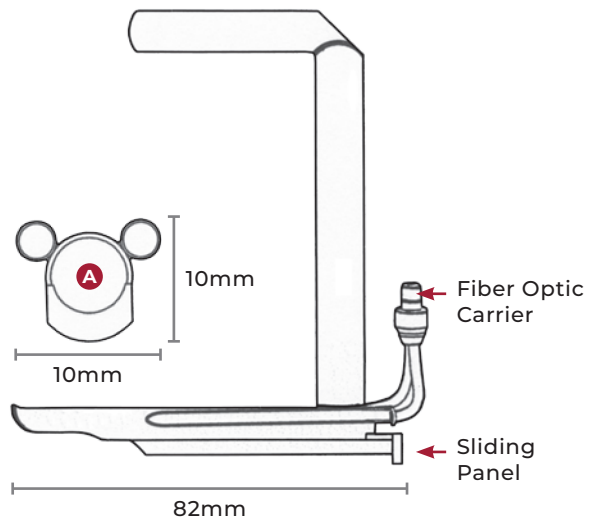

A 35370.10
NEGUS
Infant size
ø 10 x 10mm - 82mm

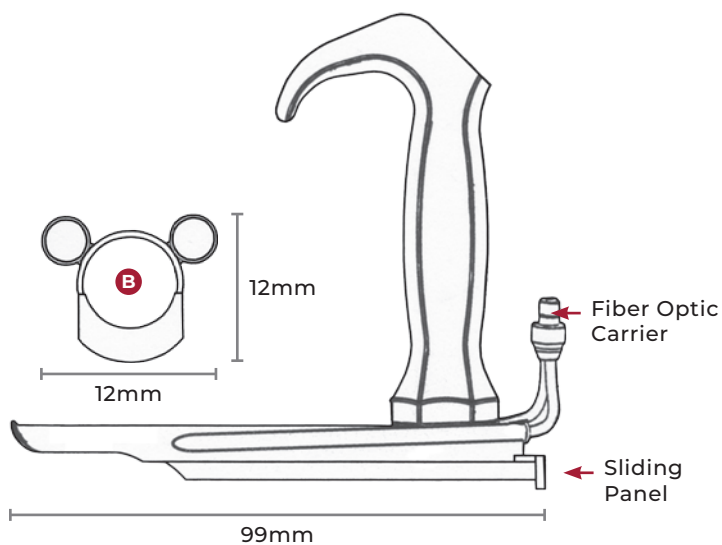


B 35370.12
NEGUS
Child size
ø 12 x 12mm - 99mm

NEGUS Laryngoscopes


Laryngoscopes with twin fiber optic light carrier and removable sliding panel. Made of polished brass.

laryngoscopes

Instruments @ 50% · Inserts to Scale

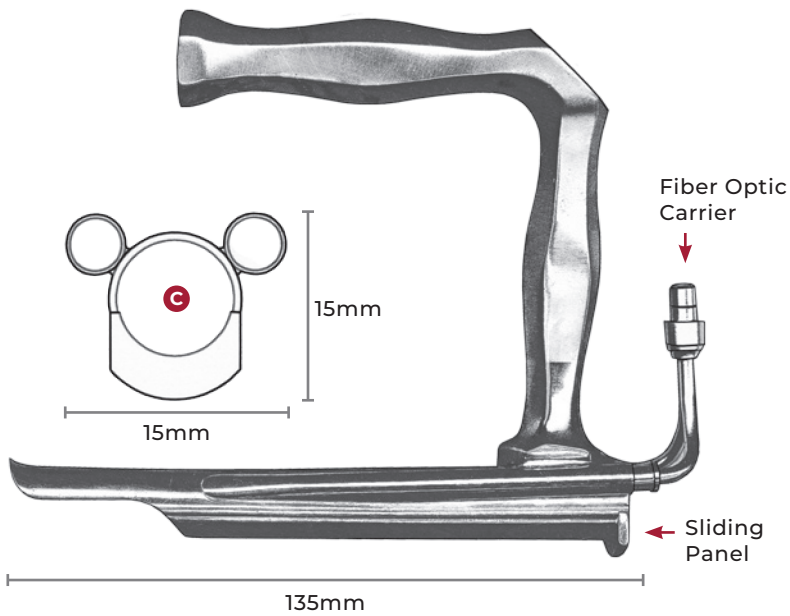

C 35370.15
NEGUS
Adult size
Ø 15 x 15mm - 135mm

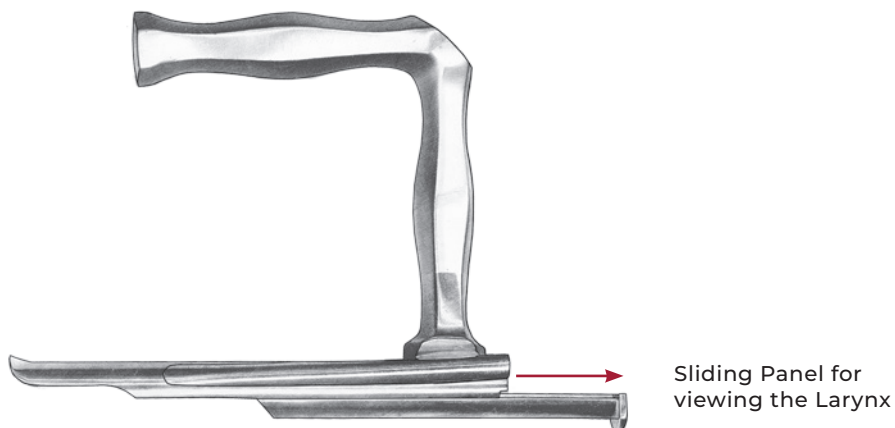

NEGUS Laryngoscopes

Laryngoscopes with twin fiber optic light carrier and removable sliding panel. Made of polished brass.

laryngoscopes


Instruments @ 50% · Inserts to Scale · * refers to instrument the displayed

Operating Laryngoscopes for 35390.01 to 35390.07



Injection Cannula
for 35390.01 to 35390.07,
fits size 1-7, 11cm

35396.11

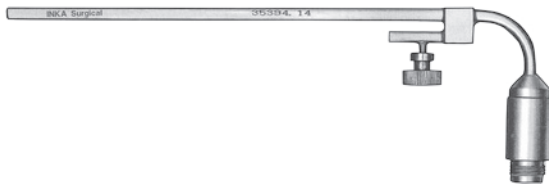

Fiber Optic Light Carriers only
for operating Laryngoscopes
35390.01 to 35390.07

35394.14* Fits size 1-2, 14cm

35394.16 Fits size 3-6, 16cm

35394.21 Fits size 7, 21cm



Vapor Suction Tube only
fits size 3-7, 15cm

35398.15

Modified for CO₂ Laser Surgery for 35392.01 to 35392.07



Fiber Optic Light Carriers only

35395.14 Fits size 1-2

35395.16* Fits size 3-6

35395.21 Fits size 7



Vapor Suction Tube
fits size 3-7, 15cm

35399.15



Twin Pattern
Fiber Optic Light Carrier only
for 35370.15

35373.10 Infant Size

35373.12 Child Size

35373.15* Adult Size




Fiber Optic Clamp
Proximal illumination with triple
connector, 11cm


16710.00

cylindrical shafts

Instruments @ 50% · Inserts to Scale

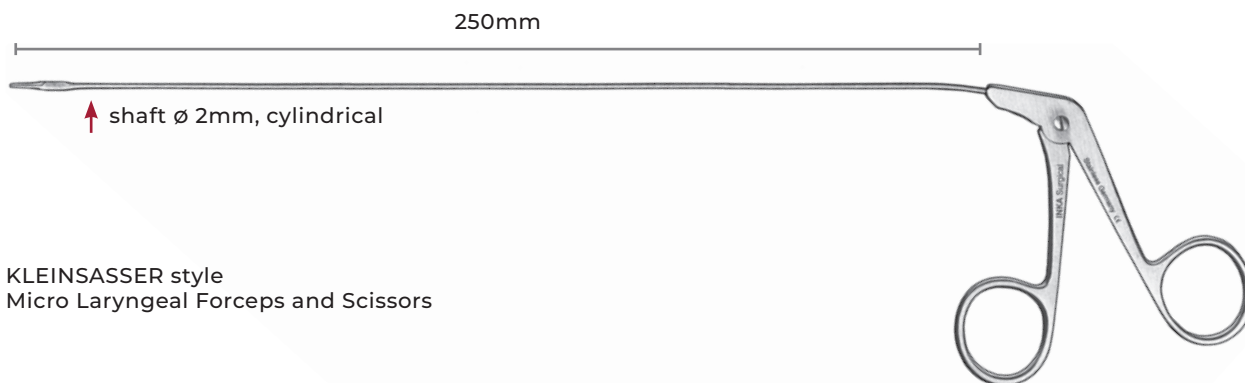


KLEINSASSER style
Micro Laryngeal Forceps and Scissors

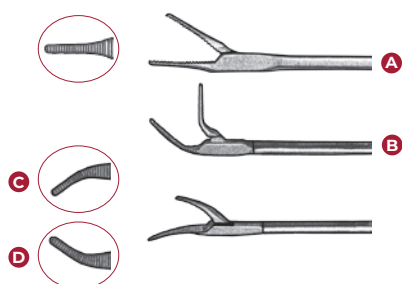


Alligator Grasping Forceps
with cylindrical shaft

- A 35430.720** Straight
- B 35430.722** Up-Curved
- C 35430.726** Left Curved
- D 35430.724** Right Curved

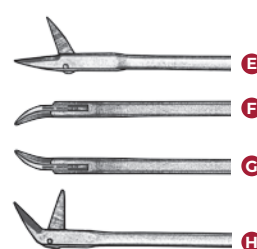


Micro Laryngeal Scissors
with cylindrical shaft

- E 35430.210** Straight
- F 35430.230** Left Curved
- G 35430.225** Right Curved
- H 35430.215** Up-Curved

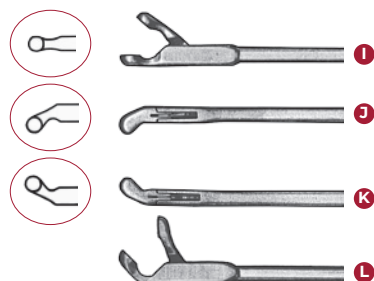


Round Cup Forceps ø 2.5mm
with cylindrical shaft

- I 35430.780** Straight
- J 35430.786** Left Curved
- K 35430.784** Right Curved
- L 35430.782** Up-Curved

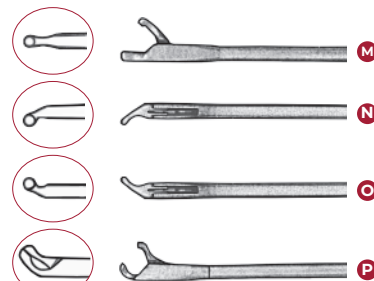


Round Cup Forceps
extra delicate cup ø 1mm
with cylindrical shaft

- M 35430.780-1** Straight
- N 35430.786-1** Left Curved
- O 35430.784-1** Right Curved
- P 35430.782-1** Up-Curved

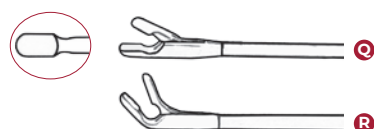

Oval Cup Forceps 3x5mm
with cylindrical shaft

- Q 35430.814** Straight
- R 35430.816** Up-Curved

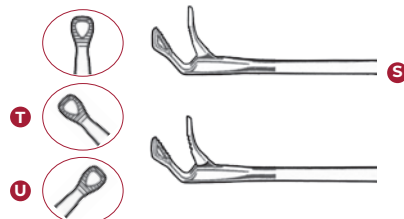


Triangular Grasping Forceps
(heart shaped) fenestrated

- S 35430.730** Up-Curved
- T 35430.734** Up-Curved 45° Left
- U 35430.733** Up-Curved 45° Right

cylindrical shafts

Instruments @ 50% · Inserts to Scale

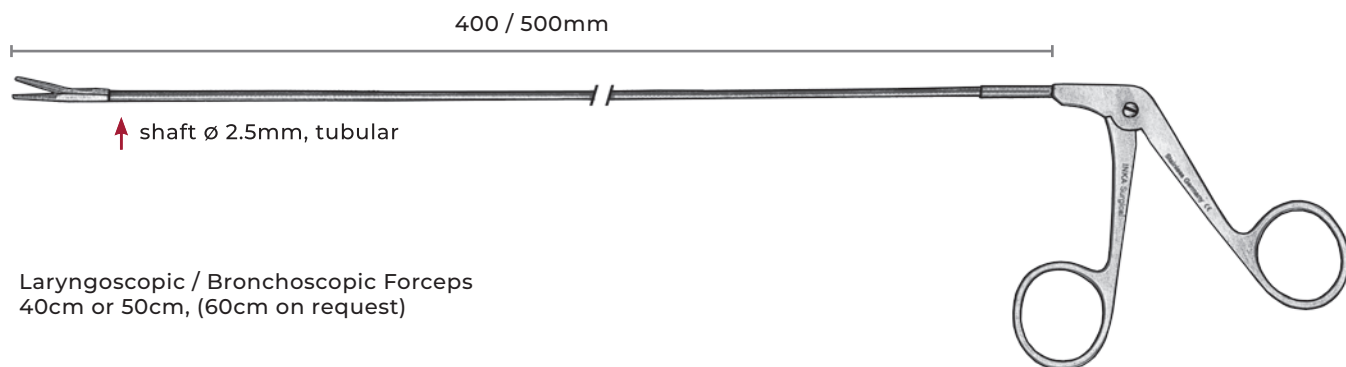

ALLIGATOR
Grasping Forceps

35271.40 40cm
35271.50 50cm



CHEVALIER JACKSON
Round Cup Forceps

35287.40 40cm
35287.50 50cm



PEANUT
Grasping Forceps

35275.40 40cm
35275.50 50cm



Vocal Nodule Forceps
up-curved

35290.40 40cm
35290.50 50cm



PIKE
Grasping Forceps

35279.40 40cm
35279.50 50cm



Grasping Forceps
with cutting tip

35294.40 40cm
35294.50 50cm



PATTERSON
Oval Cup Forceps

35283.40 40cm
35283.50 50cm




CHEVALIER JACKSON
Basket Punch

35299.40 40cm
35299.50 50cm

cylindrical shafts

Instruments @ 50% · Inserts to Scale

350 / 500mm

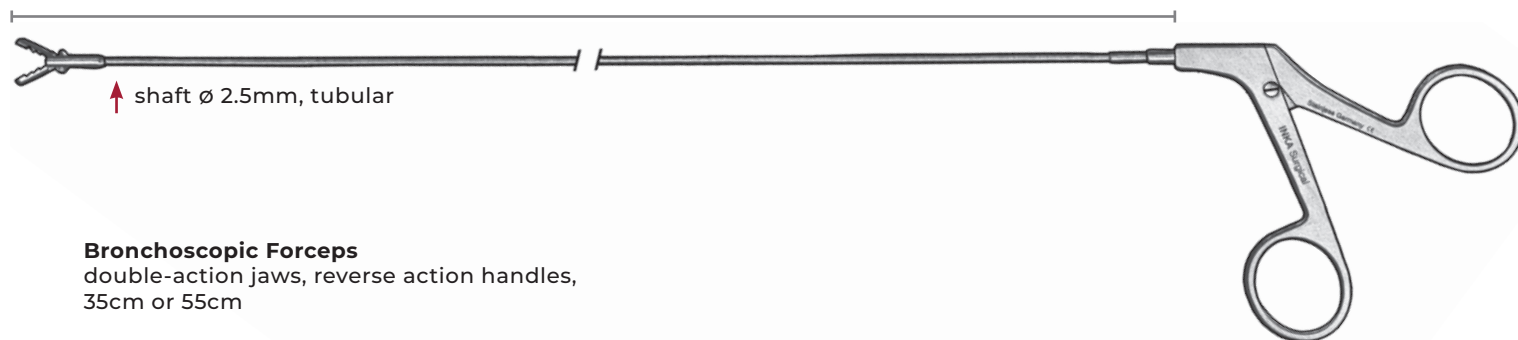

Bronchoscopic Forceps

double-action jaws, reverse action handles,
35cm or 55cm



PIKE
Bronchoscopic Forceps

35310.35 35cm
35310.55 55cm



UNIVERSAL
Foreign Body Forceps

35314.35 35cm
35314.55 55mm



Biopsy & Grasping
round cup

35318.35 35cm
35318.55 55cm





Grasping Forceps
with fenestrated jaws

35321.35 35cm
35321.55 55mm


conical shafts

Instruments @ 50% · Inserts to Scale

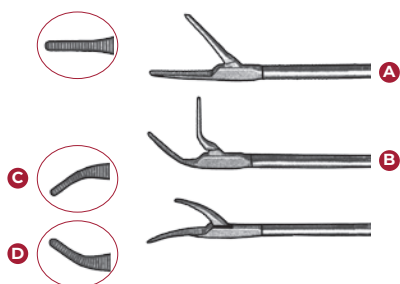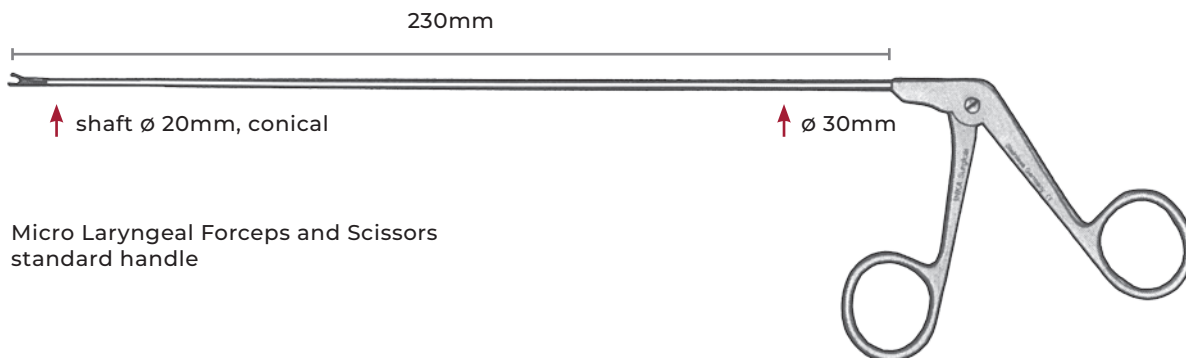


Alligator Grasping Forceps

- A 35420.720** Straight
- B 35420.722** Up-Curved
- C 35420.726** Left Curved
- D 35420.724** Right Curved

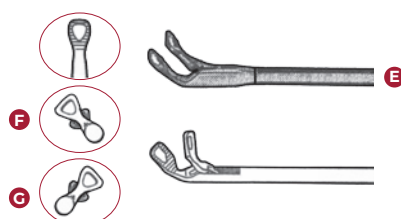


Triangular Grasping Forceps (heart shaped) fenestrated

- E 35420.730** Up-Curved
- F 35420.734** Up-Curved 45° Left
- G 35420.733** Up-Curved 45° Right

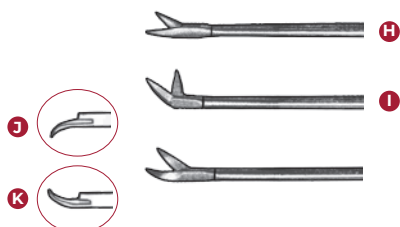


Micro Laryngeal Scissors

- H 35420.210** Straight
- I 35420.215** Up-Curved
- J 35420.230** Left Curved
- K 35420.225** Right Curved

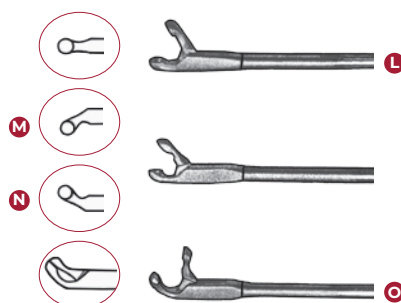

Round Cup Forceps $\varnothing$ 2.5mm cups

- L 35420.780** Straight
- M 35420.786** Left Curved
- N 35420.784** Right Curved
- O 35420.782** Up-Curved


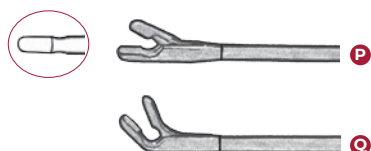


Oval Cup Forceps 2.5x4mm cup


- P 35420.814** Straight
- Q 35420.816** Up-Curved

reverse action handles

Instruments @ 50% · Inserts to Scale

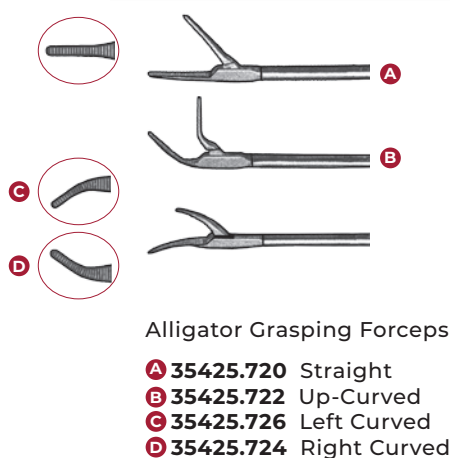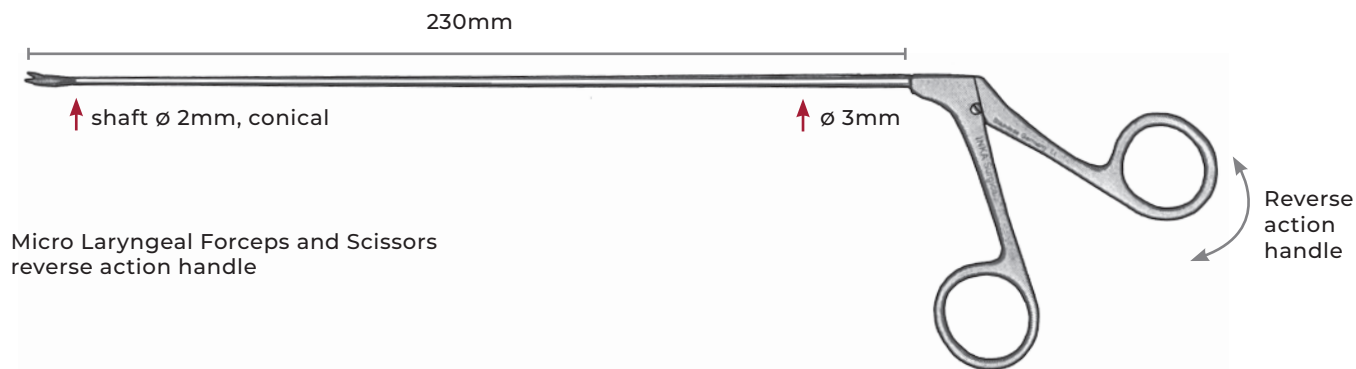


Alligator Grasping Forceps

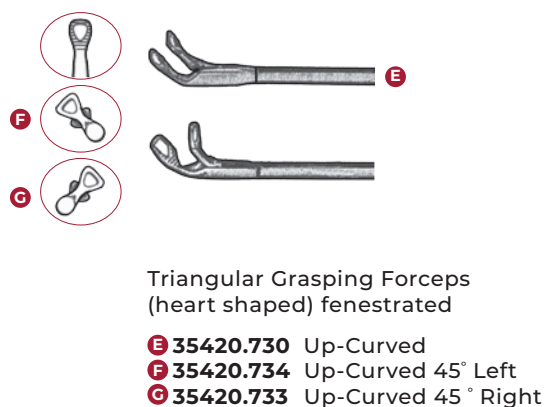


Triangular Grasping Forceps (heart shaped) fenestrated



Micro Laryngeal Scissors

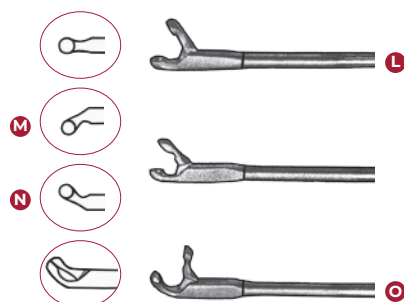

Round Cup Forceps ø 2.5mm cup

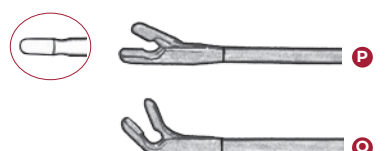

Oval Cup Forceps 2.5x4mm cup

L 35425.780 Straight

M 35425.786 Left Curved

N 35425.784 Right Curved


O 35425.782 Up-Curved

P 35425.814 Straight


Q 35425.816 Up-Curved

paediatric pattern

Instruments @ 50% · Inserts to Scale

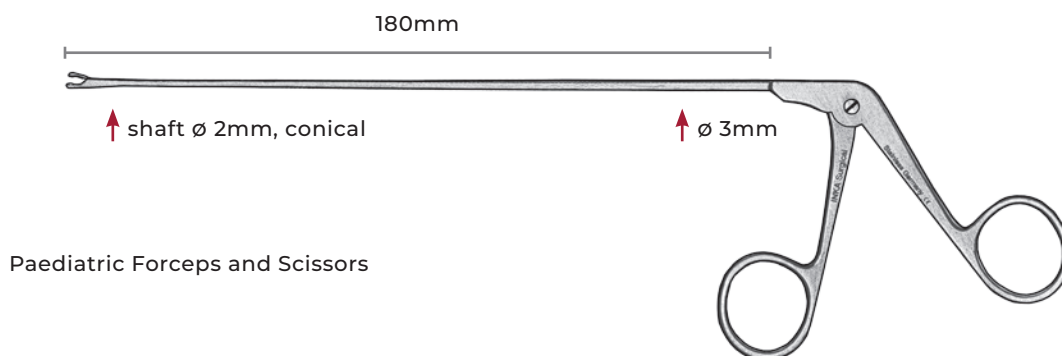

Paediatric Forceps and Scissors

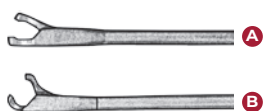


KLEINSASSER style
Paediatric Micro Cup Forceps
 $\varnothing$ 1mm cup

- A 35422.780-1** Straight
- B 35422.782-1** Up-Curved



KLEINSASSER style
Paediatric Alligator Forceps
straight

35422.720

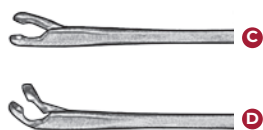


KLEINSASSER style
Paediatric Cup Forceps
 $\varnothing$ 2mm cup

- C 35422.780-2** Straight
- D 35422.782-2** Up-Curved



KLEINSASSER style
Paediatric Micro Laryngeal
Scissors

- G 35422.210** Straight
- H 35422.230** Left Curved
- I 35422.225** Right Curved
- J 35422.215** Up-Curved

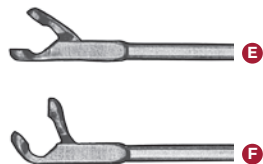

KLEINSASSER style
Paediatric Cup Forceps
 $\varnothing$ 3mm cup

- E 35422.780-3** Straight
- F 35422.782-3** Up-Curved




Round Cup Forceps
with cylindrical shaft $\varnothing$ 2.5mm,
 $\varnothing$ 4mm cup, straight, 18cm

33928.18

paediatric pattern

Instruments @ 50% · Inserts to Scale

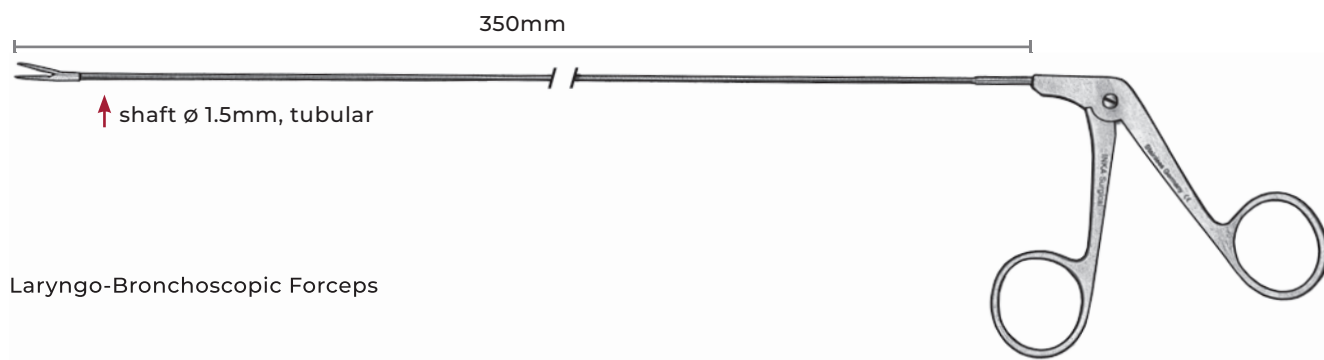

ALLIGATOR
Grasping Forceps
35240.35



CHEVALIER JACKSON
Round Cup Forceps
35256.35



PEANUT
Grasping Forceps
35244.35



Vocal Nodule Forceps
up-curved
35259.35



PIKE
Grasping Forceps
35248.35




Grasping Forceps
with cutting tip
35263.35



PATERSON
Oval Cup Forceps
35252.35

paediatric pattern

Instruments @ 50% · Inserts to Scale

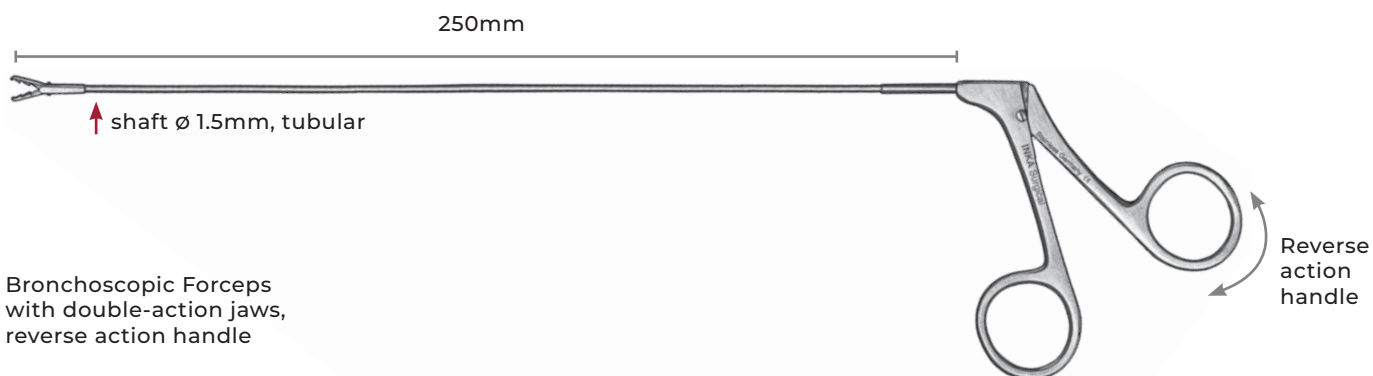

PIKE
Foreign Body Forceps
35302.25



Biopsy and Grasping Forceps
35304.25




Biopsy Round Cup Forceps
round cup $\varnothing$ 2mm
35306.25




Grasping Forceps
with fenestrated jaws
35308.25

miniature range

Instruments @ 50% · Inserts to Scale

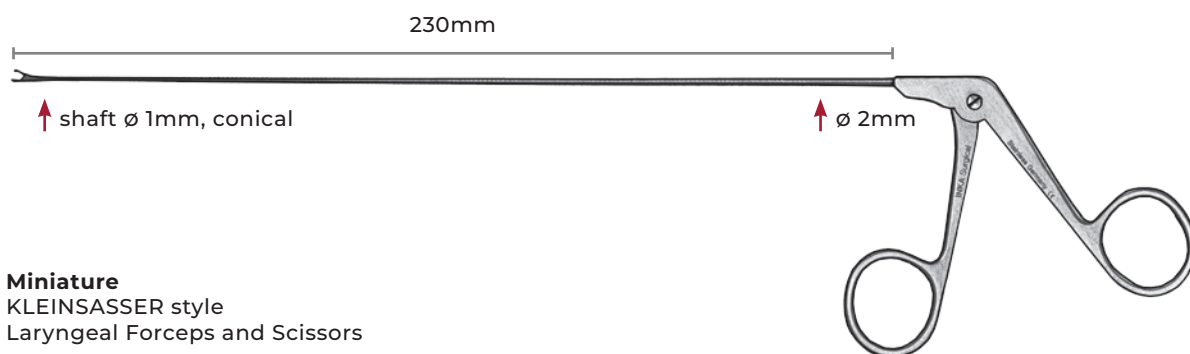


Miniature
KLEINSASSER style
Laryngeal Forceps and Scissors

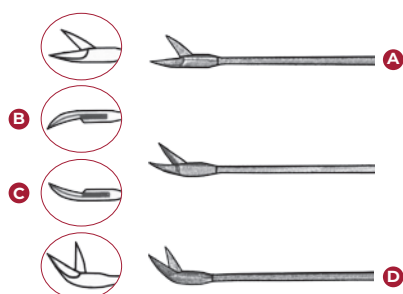


Miniature Laryngeal Scissors

- A 35450.212** Straight
- B 35450.232** Left Curved
- C 35450.227** Right Curved
- D 35450.217** Up-Curved

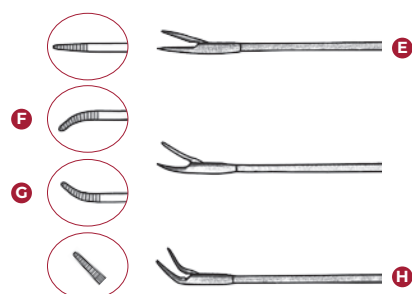


Miniature Alligator
Grasping Forceps

- E 35450.910** Straight
- F 35450.916** Left Curved
- G 35450.914** Right Curved
- H 35450.912** Up-Curved

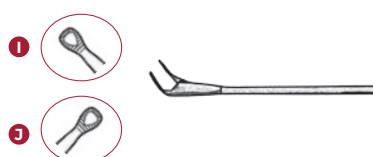

Miniature Grasping Forceps
with triangular jaws,
(heart shaped) fenestrated,
45° up-curved,

- I 35450.928** Right
- J 35450.929** Left

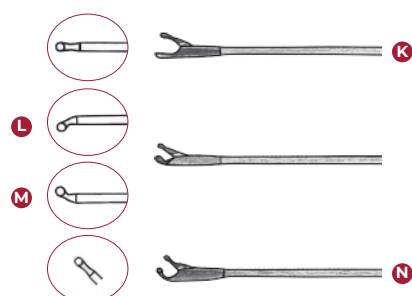


Miniature Round Cup Forceps
ø 1mm cup


- K 35450.930** Straight
- L 35450.936** Left Curved
- M 35450.934** Right Curved
- N 35450.932** Up-Curved

miniature range

Instruments @ 50% · Inserts to Scale

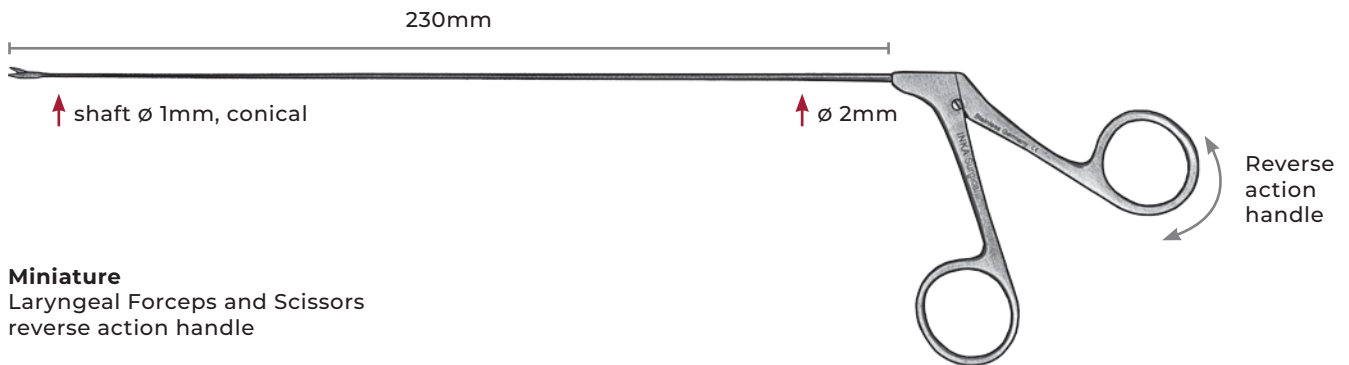


Miniature
Laryngeal Forceps and Scissors
reverse action handle

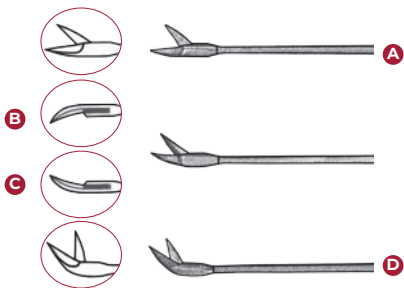


Miniature Laryngeal Scissors

- A 35455.212** Straight
- B 35455.232** Left Curved
- C 35455.227** Right Curved
- D 35455.217** Up-Curved

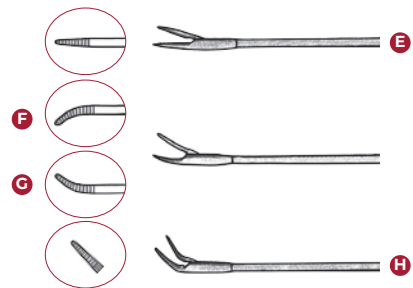


Miniature Alligator
Grasping Forceps

- E 35455.910** Straight
- F 35455.916** Left Curved
- G 35455.914** Right Curved
- H 35455.912** Up-Curved

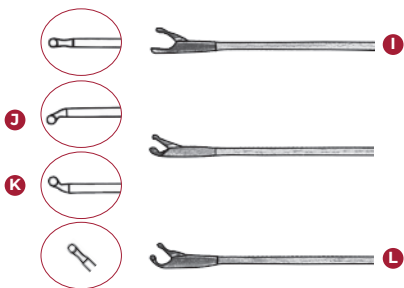

Miniature Round Cup Forceps
Ø 1mm cup

- I 35455.930** Straight
- J 35455.936** Left Curved
- K 35455.934** Right Curved
- L 35455.932** Up-Curved

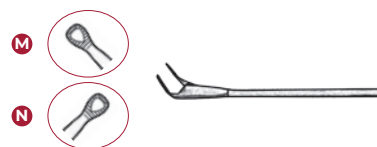


Miniature Grasping Forceps
with triangular jaws,
(heart shaped) fenestrated


- M 35455.929** Up-Curved 45° Left
- N 35455.928** Up-Curved 45° Right

cleaning channel range

Instruments @ 50% · Inserts to Scale

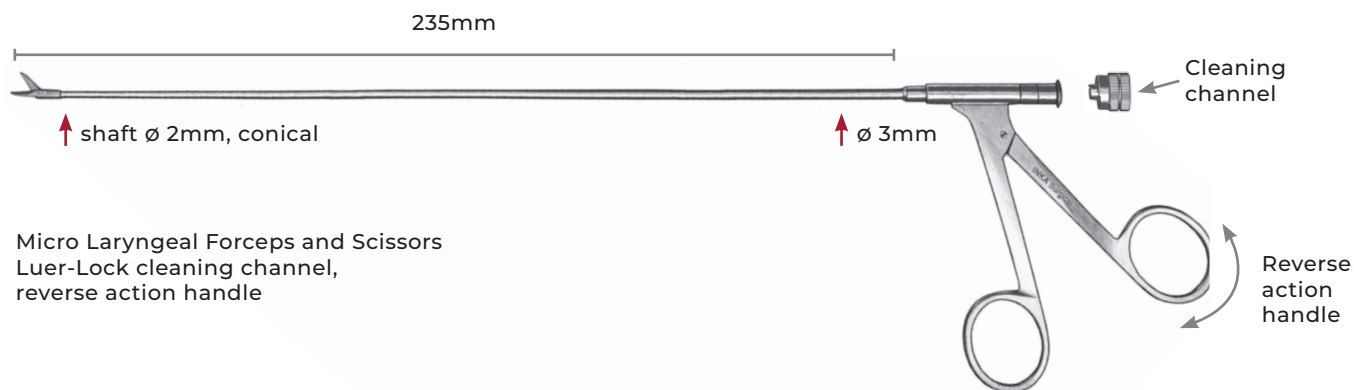


Micro Laryngeal Forceps and Scissors
Luer-Lock cleaning channel,
reverse action handle

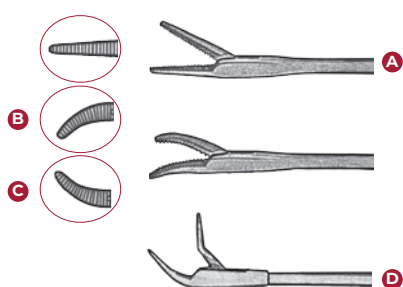


Alligator Grasping Forceps

- A 35461.610** Straight
- B 35461.612** Left Curved
- C 35461.613** Right Curved
- D 35461.611** Up-Curved

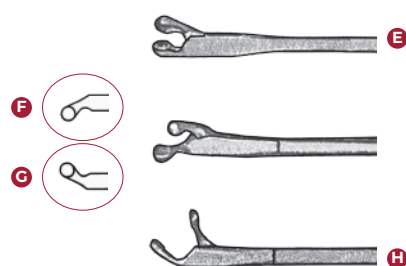


Round Cup Forceps
 $\varnothing$ 2mm cup

- E 35461.600** Straight
- F 35461.602** Left Curved
- G 35461.603** Right Curved
- H 35461.601** Up-Curved

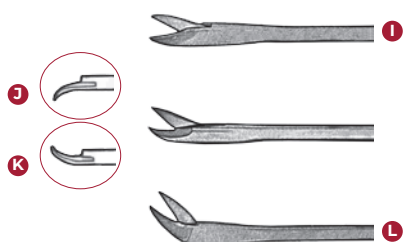


Laryngeal Micro Scissors

- I 35461.621** Straight
- J 35461.623** Left Curved
- K 35461.624** Right Curved
- L 35461.622** Up-Curved

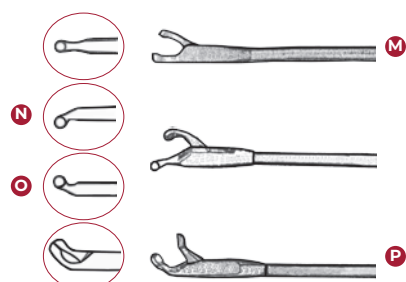

Round Cup Forceps
 $\varnothing$ 1mm cup

- M 35461.640** Straight
- N 35461.642** Left Curved
- O 35461.643** Right Curved
- P 35461.641** Up-Curved

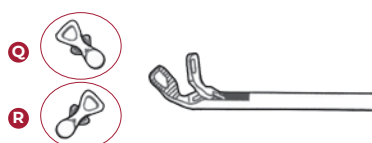

Laryngeal Grasping Forceps
with triangular jaws,
(heart shaped) fenestrated,
45° up-curved

- Q 35461.630** Left
- R 35461.631** Right


cleaning channel range

Instruments @ 50% · Inserts to Scale

270mm




Micro Laryngeal Forceps and Scissors
Luer-Lock cleaning channel,
reverse action handle

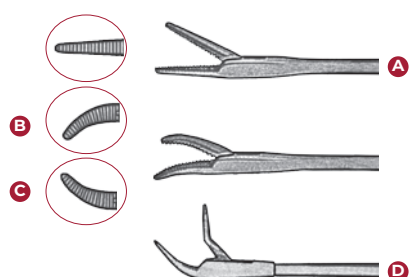


Alligator Grasping Forceps

- A 35464.670** Straight
- B 35464.672** Left Curved
- C 35464.673** Right Curved
- D 35464.671** Up-Curved

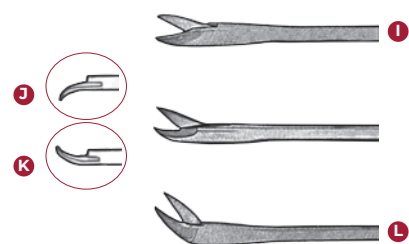


Laryngeal Micro Scissors

- I 35464.681** Straight
- J 35464.683** Left Curved
- K 35464.684** Right Curved
- L 35464.682** Up-Curved

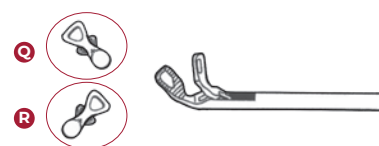


Laryngeal Grasping Forceps
with triangular jaws,
(heart shaped) fenestrated,
45° up-curved,

- Q 35464.690** Left
- R 35464.691** Right

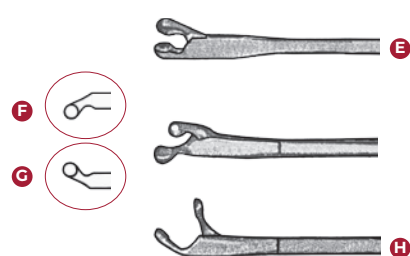

Round Cup Forceps
ø 2mm cup

- E 35464.660** Straight
- F 35464.662** Left Curved
- G 35464.663** Right Curved
- H 35464.661** Up-Curved

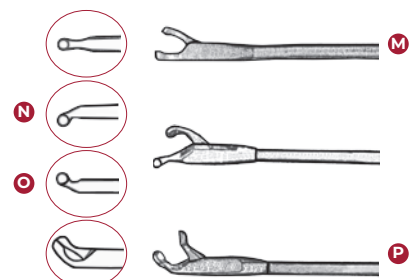


Round Cup Forceps
ø 1mm cup

- M 35464.650** Straight
- N 35464.652** Left Curved
- O 35464.653** Right Curved
- P 35464.651** Up-Curved


cup forceps - scissors

Instruments @ 50% · Inserts to Scale

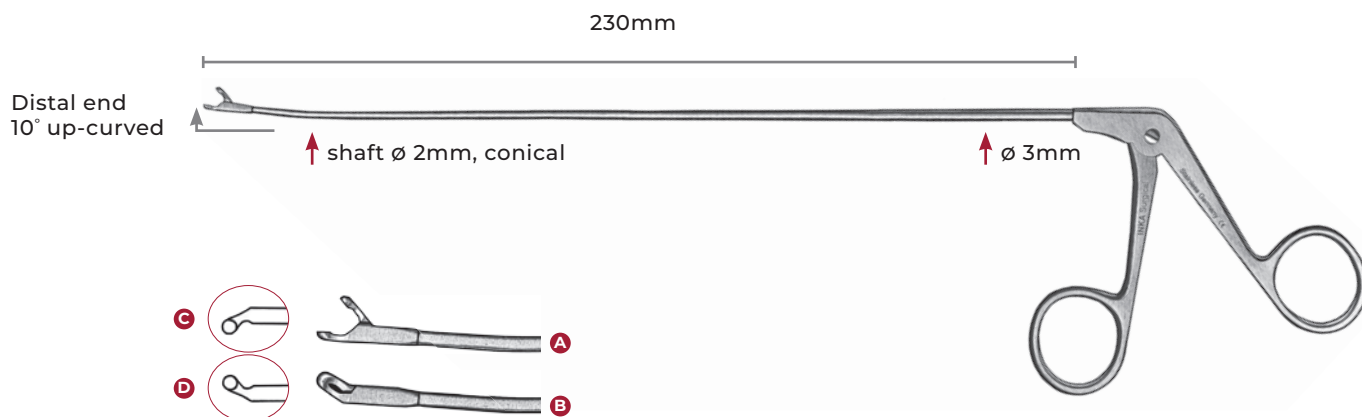

KLEINSASSER style
Micro Laryngeal Cup Forceps
 $\varnothing$ 2mm cup

- A** **35470.1107** Straight
- B** **35470.1108** Up-Curved
- C** **35470.1109** Right Curved
- D** **35470.1110** Left Curved

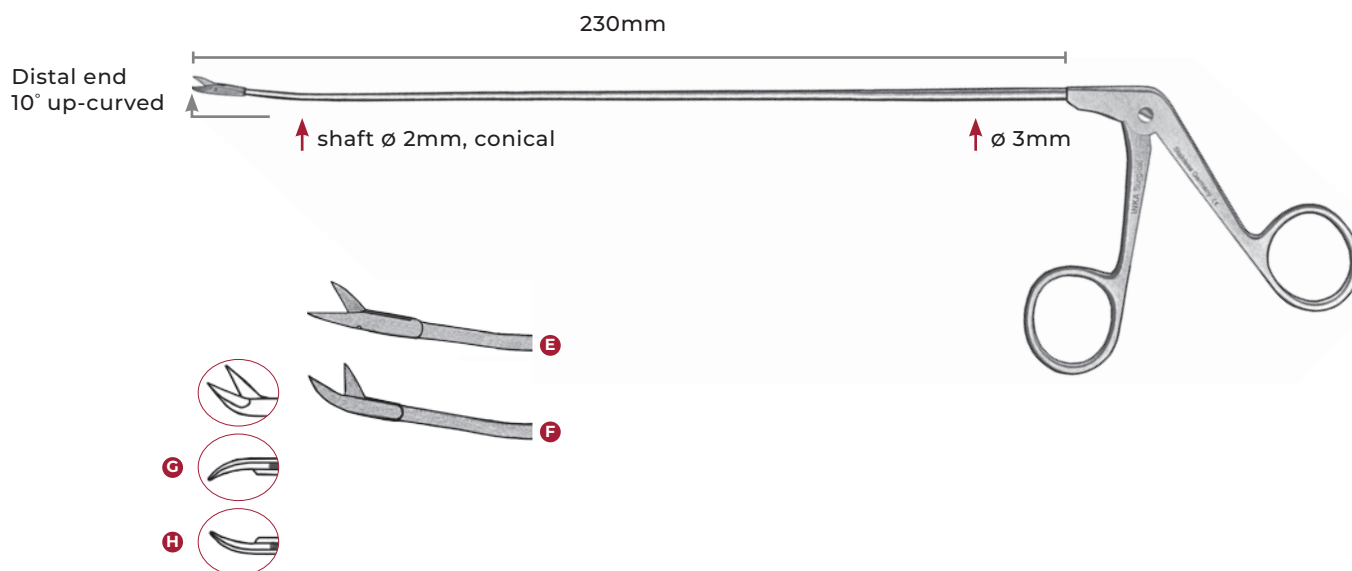


KLEINSASSER style
Micro Laryngeal Scissors

- E** **35472.1111** Straight
- F** **35472.1112** Up-Curved
- G** **35472.1114** Left Curved
- H** **35472.1113** Right Curved

cup forceps - scissors

Instruments @ 50% · Inserts to Scale

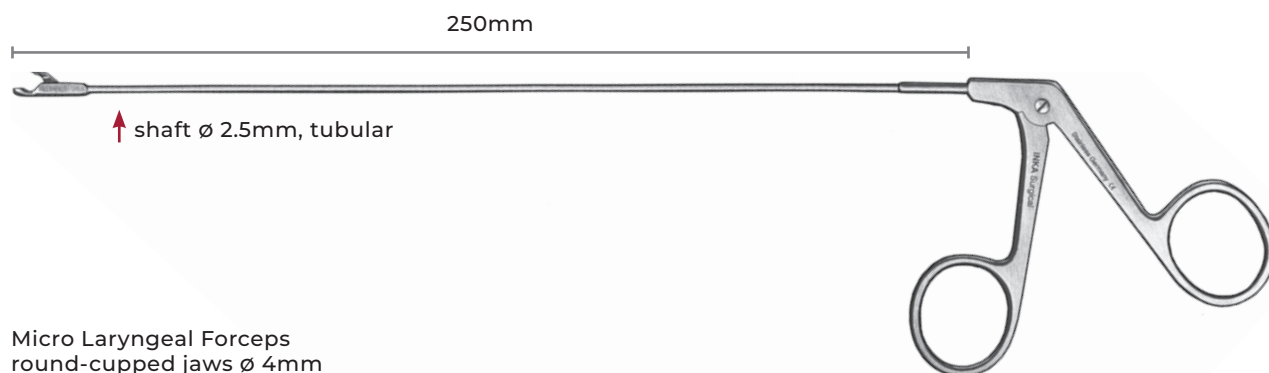

Round Cup Forceps
ø 4mm cup, straight
35480.04



Round Cup Forceps
ø 4mm cup, right curved
35483.04





Round Cup Forceps
ø 4mm cup, left curved
35482.04



Round Cup Forceps
ø 4mm cup, up-curved
35486.04



needleholder - vocal cord spreader

Instruments @ 50% · Inserts to Scale

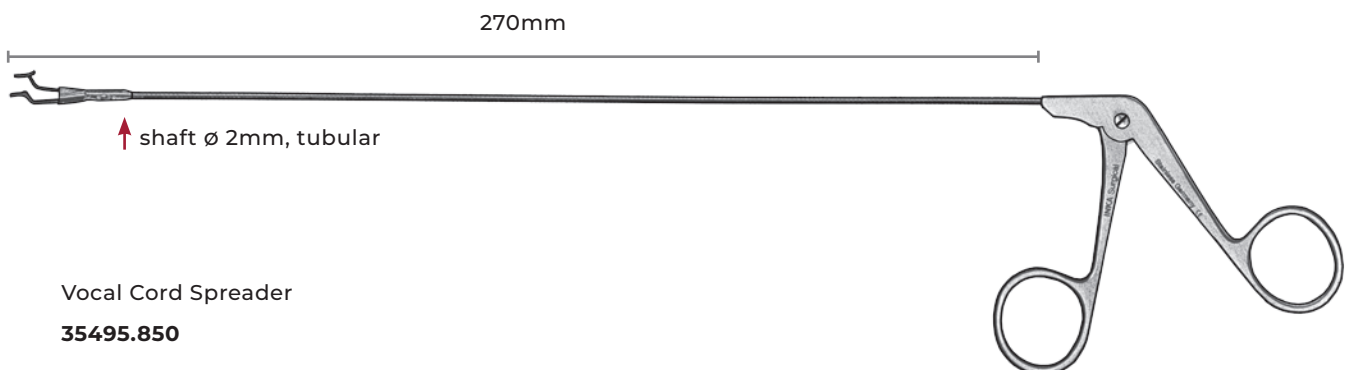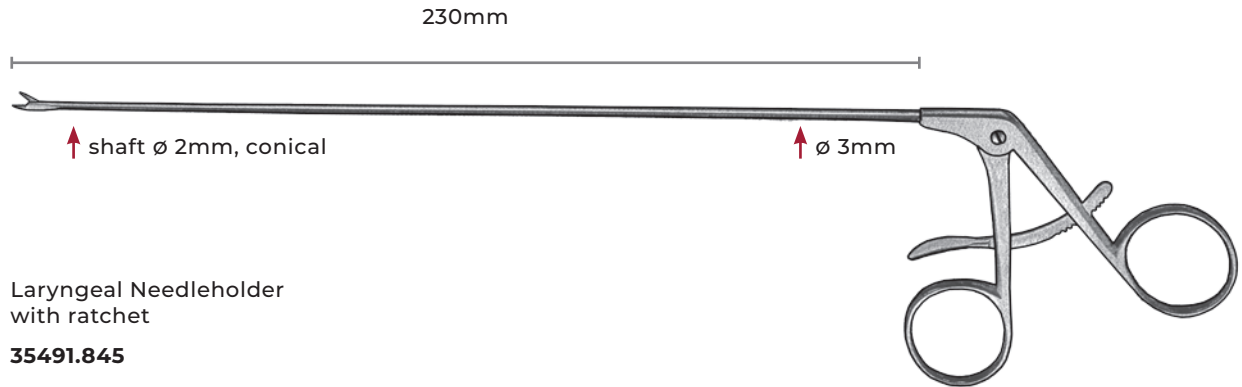


foreign body forceps

Instruments @ 50% · Inserts to Scale

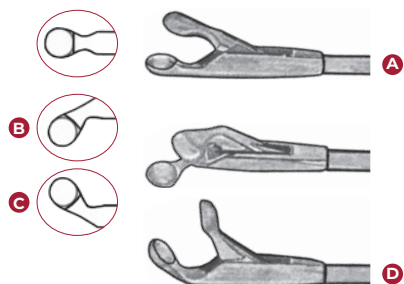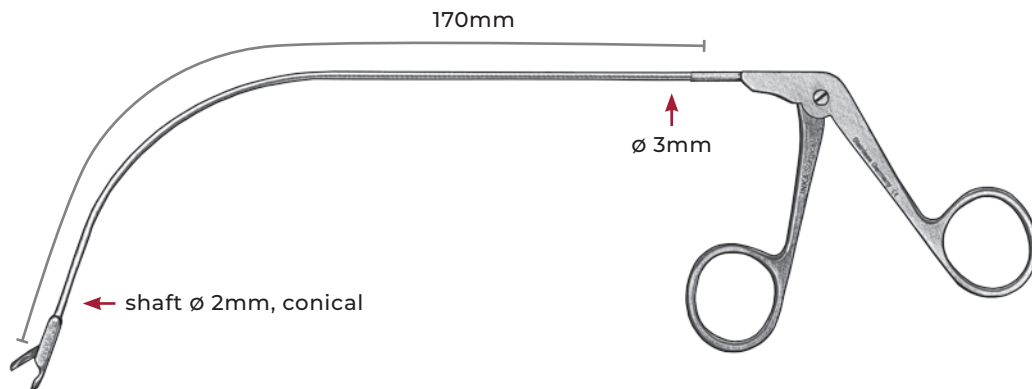

Laryngeal Biopsy Forceps
round tips ø 4mm

- A 35150.876** Straight jaw
- B 35150.884** Left Curved
- C 35150.882** Right Curved
- D 35150.878** Up-Curved

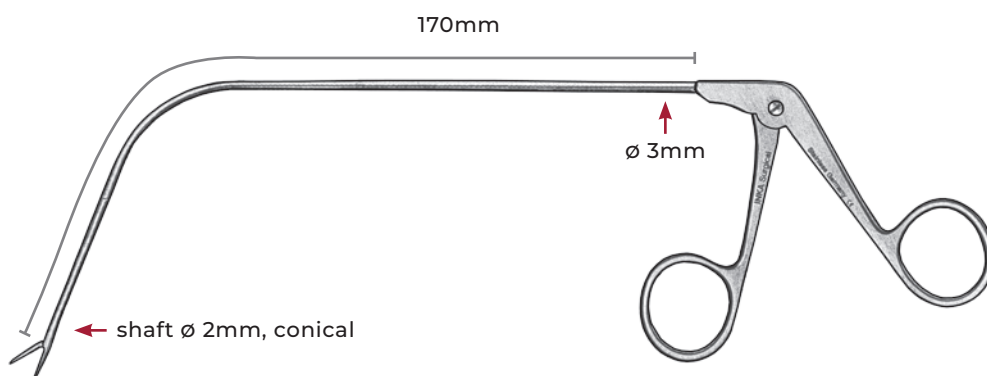


Laryngeal Alligator Forceps
serrated, for foreign body
removal

35153.862

foreign body forceps

Instruments @ 50% · Inserts to Scale

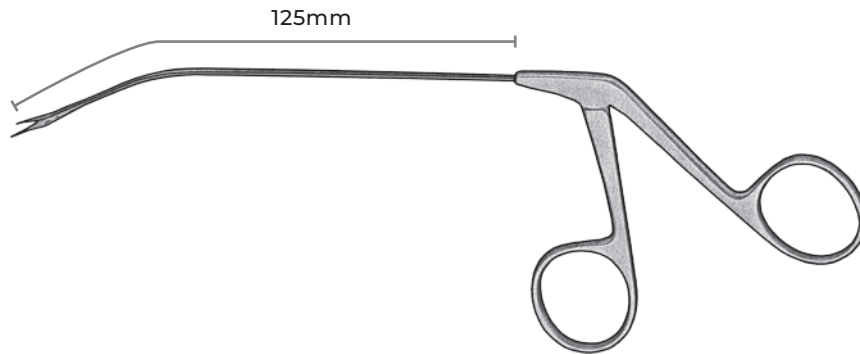


KALUSKAR
Fish Bone Removing Forceps
opening to right

35156.12

esophagoscopy forceps

Instruments @ 50% · Inserts to Scale

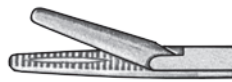

ALLIGATOR
Grasping Forceps

35174.30 30cm
35174.40 40cm
35174.50 50cm



PATERSON
Oval Cup Forceps

35186.30 30cm
35186.40 40cm
35186.50 50cm




PIKE
Grasping Forceps

35178.30 30cm
35178.40 40cm
35178.50 50cm



JACKSON
Round Cup Forceps
ø 4mm cup

35190.30 30cm
35190.40 40cm
35190.50 50cm

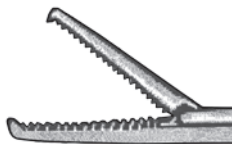

JACKSON
Rotation Forceps
with 1x2 teeth

35182.30 30cm
35182.40 40cm
35182.50 50cm




Vocal Nodule Forceps
up-curved, ø 4mm cup

35194.30 30cm
35194.40 40cm
35194.50 50cm

esophagoscopy forceps

Instruments @ 50% · Inserts to Scale

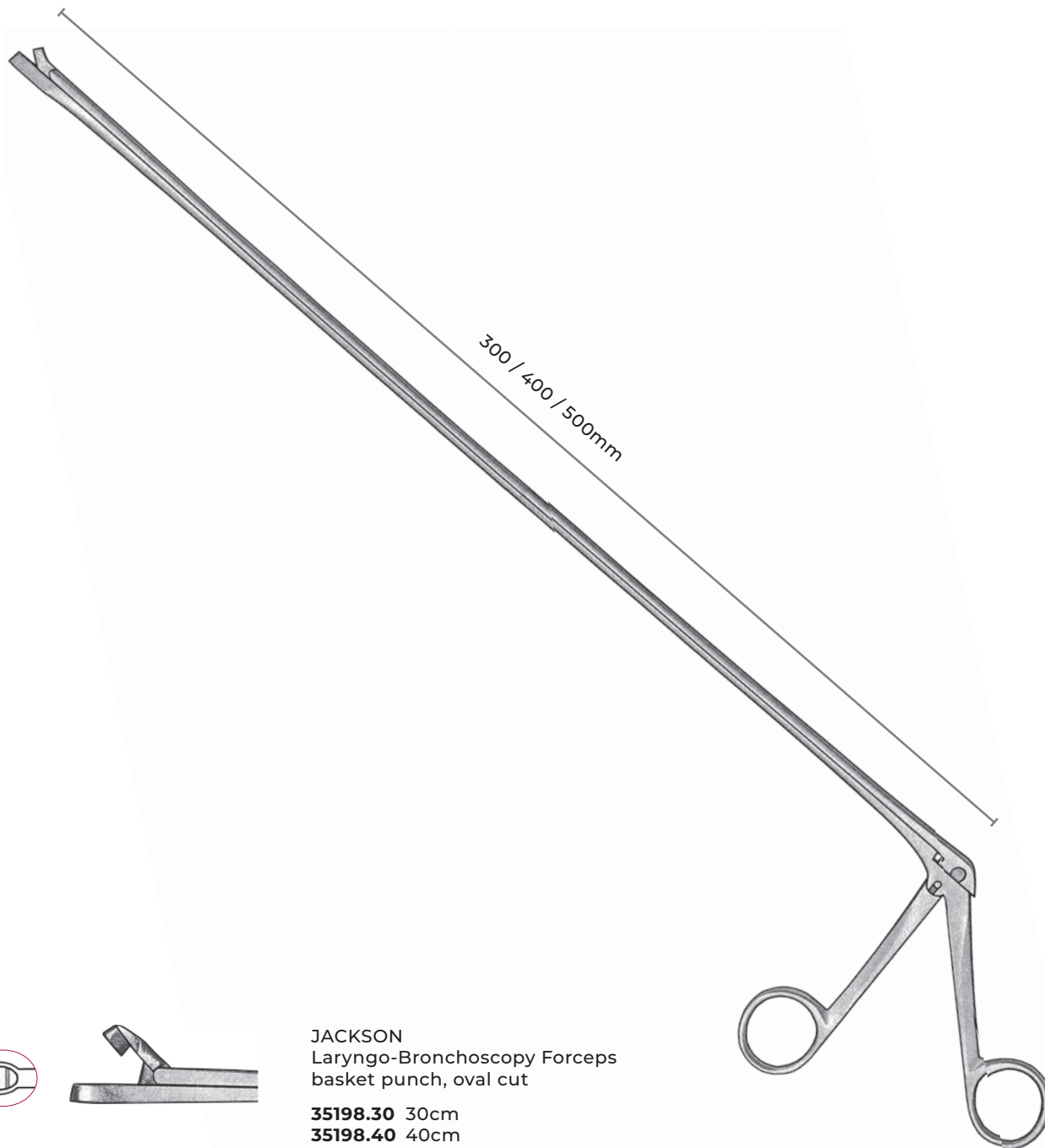


JACKSON
Laryngo-Bronchoscopy Forceps
basket punch, oval cut

35198.30 30cm

35198.40 40cm

35198.50 50cm

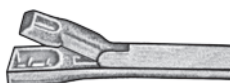

JACKSON
Laryngo-Bronchoscopy Forceps
basket punch, square cut












35202.30 30cm

35202.40 40cm

35202.50 50cm




hand instruments

Instruments @ 50% · Inserts to Scale

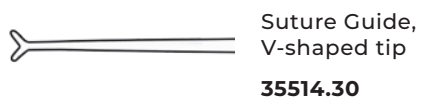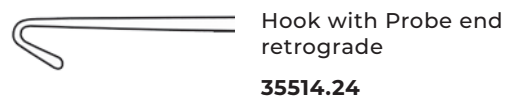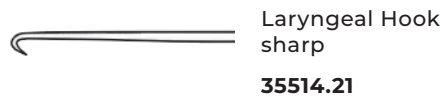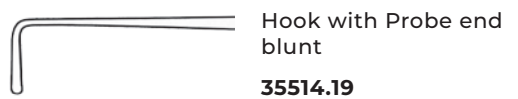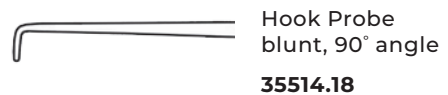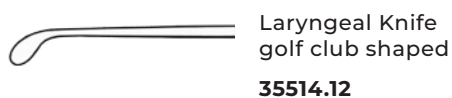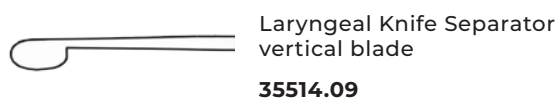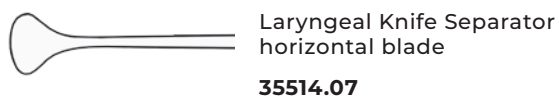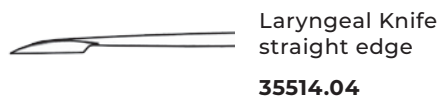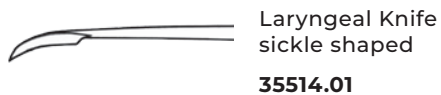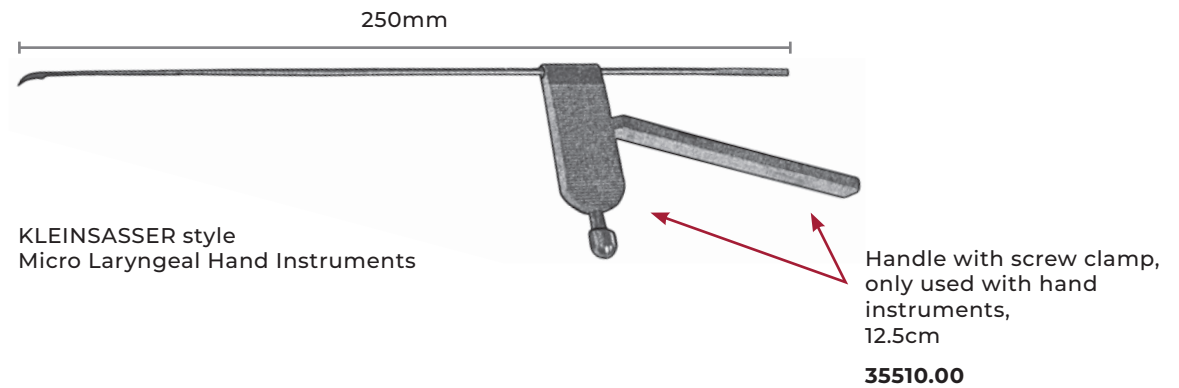

suction tubes

Instruments @ 50% · Inserts to Scale

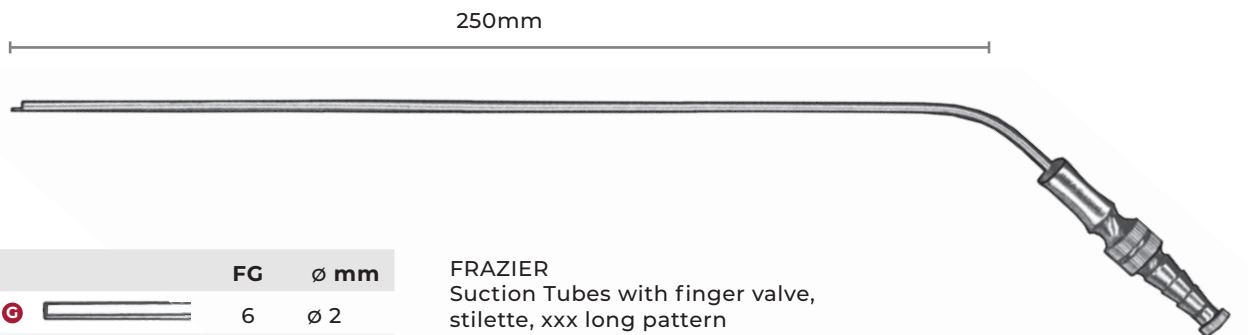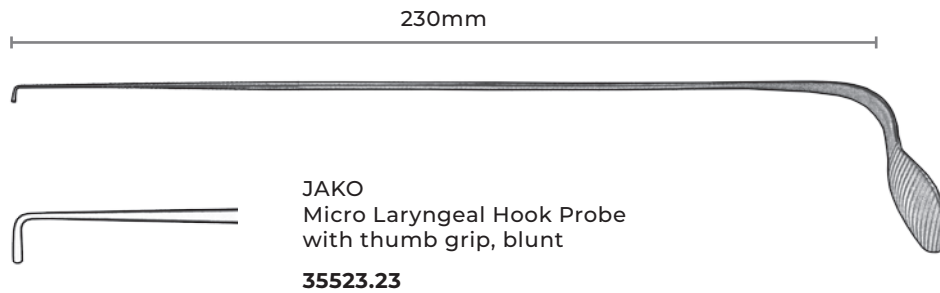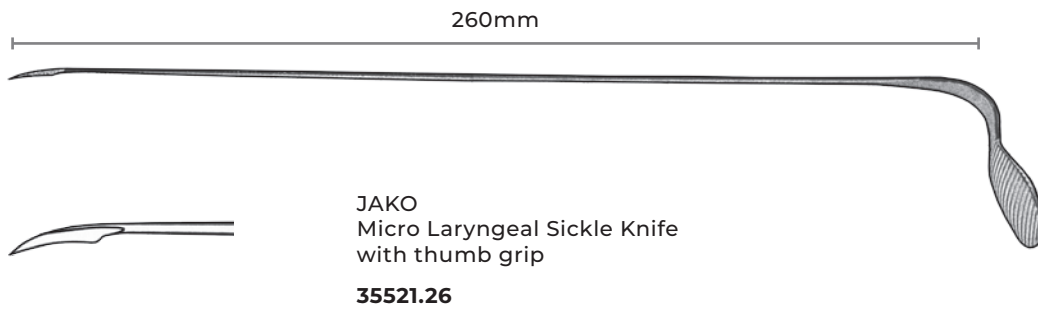

	FG	ø mm
G	6	ø 2
	7	ø 2.3
H	8	ø 2.7
I	9	ø 3
J	10	ø 3.3
K	12	ø 4
L	15	ø 5

- G 15709.06** ø 2mm
- H 15709.08** ø 2.7mm
- I 15709.09** ø 3mm
- J 15709.10** ø 3.3mm
- K 15709.12** ø 4mm
- L 15709.15** ø 5mm

15709.99 Spare Stilettes only
(nominate FG size)

mirrors

Instruments @ 50% · Inserts to Scale



Size 0 -
ø 10mm



Size 1 -
ø 12mm




Size 2 -
ø 14mm



Size 3 -
ø 16mm




Size 4 -
ø 18mm

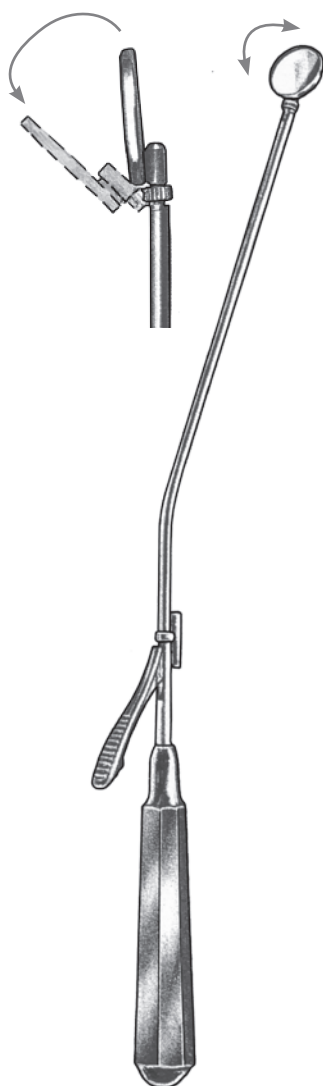

MICHEL
Rhinoscopic Mirrors
not interchangeable,
adjustable
28cm

33782.01 Size 1
33782.02 Size 2
33782.03 Size 3
33782.04 Size 4


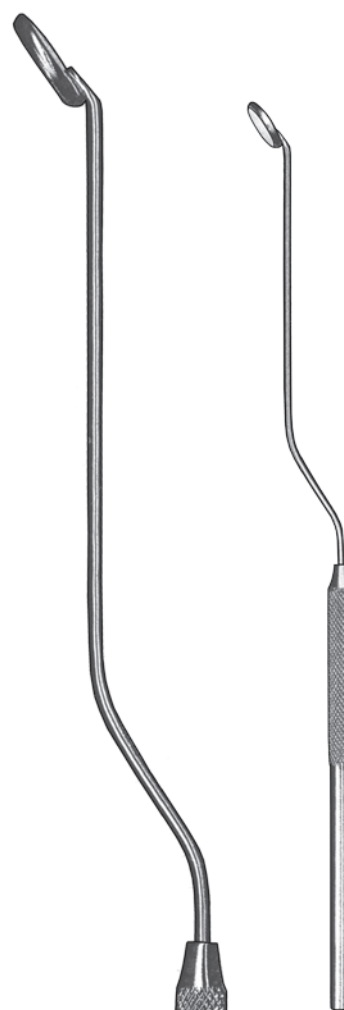


ST. CLAIR THOMSON
Laryngeal Mirrors,
cranked shaft
24cm

33778.10 Size 0
33778.12 Size 1
33778.14 Size 2
33778.16 Size 3

mirrors

Instruments @ 50% · Inserts to Scale

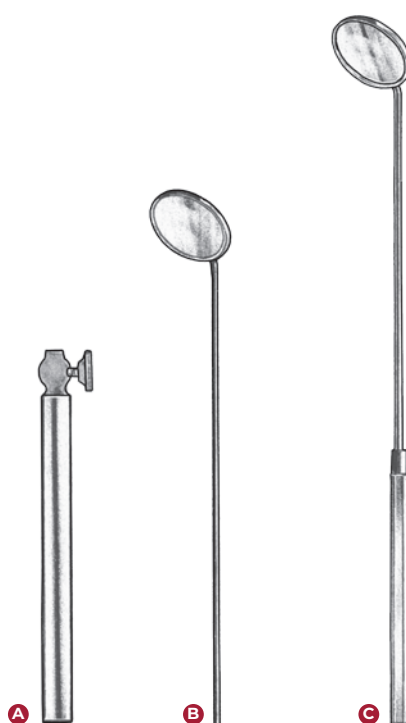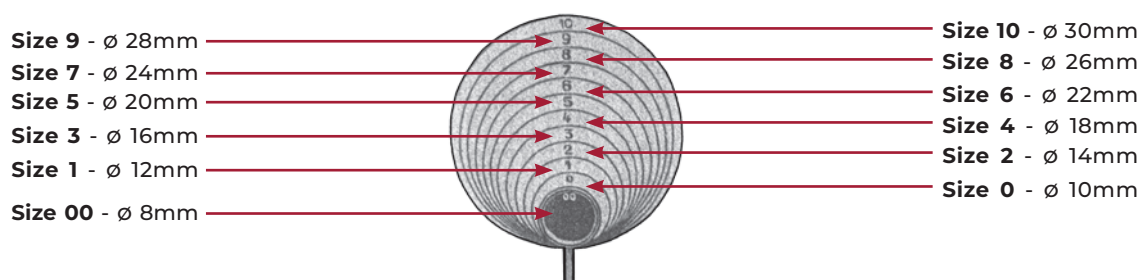

A 35104.00
Handle for mirrors
with stem chrome plated,
10cm

B Laryngeal Mirrors
with stem,
13cm

35101.08 Size 00
35101.10 Size 0
35101.12 Size 1
35101.14 Size 2
35101.16 Size 3
35101.18 Size 4
35101.20 Size 5
35101.22 Size 6
35101.24 Size 7
35101.26 Size 8
35101.28 Size 9
35101.30 Size 10

C Laryngeal Mirrors
with hexagonal handle,
18cm


35109.08 Size 00
35109.10 Size 0
35109.12 Size 1
35109.14 Size 2
35109.16 Size 3
35109.18 Size 4
35109.20 Size 5
35109.22 Size 6
35109.24 Size 7
35109.26 Size 8
35109.28 Size 9
35109.30 Size 10

laryngeal applicators

Instruments @ 50% · Inserts to Scale




Laryngeal
Cotton Applicator
straight,
28cm
35114.28

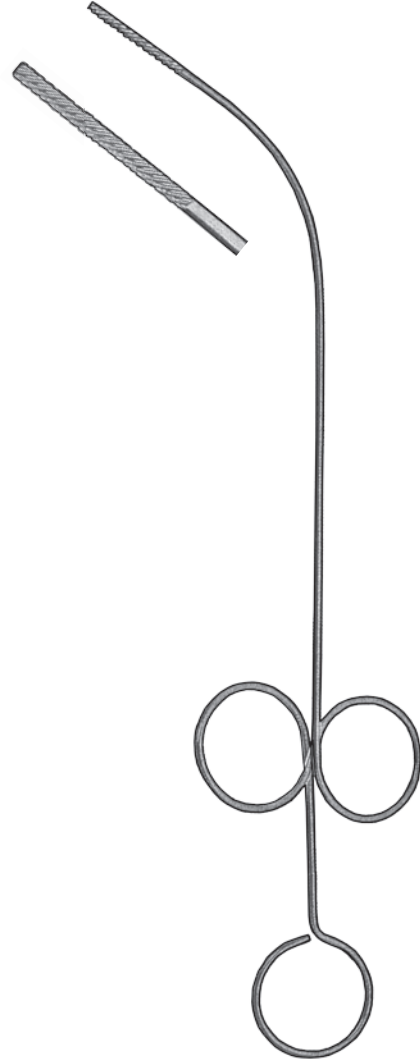

FEIN
Laryngeal Cotton Applicator
28cm
35118.28


scissors - forceps - retractors

Instruments @ 50% · Inserts to Scale

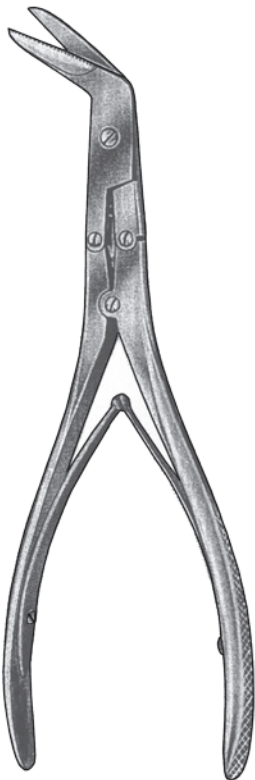

KILLIAN
Laryngectomy Scissors
angled to side,
21cm

35132.21

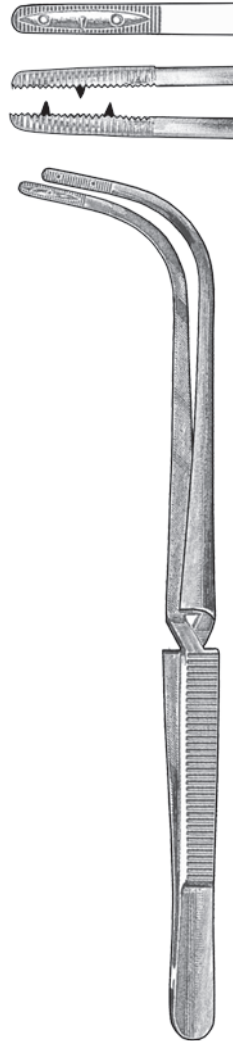

JACKSON
Laryngeal Forceps
cross-action with pins,
24cm

35136.24




REICHERT
Epiglottis Retractor
24cm

35140.24

SUTTER monopolar forceps

Instruments are not to scale

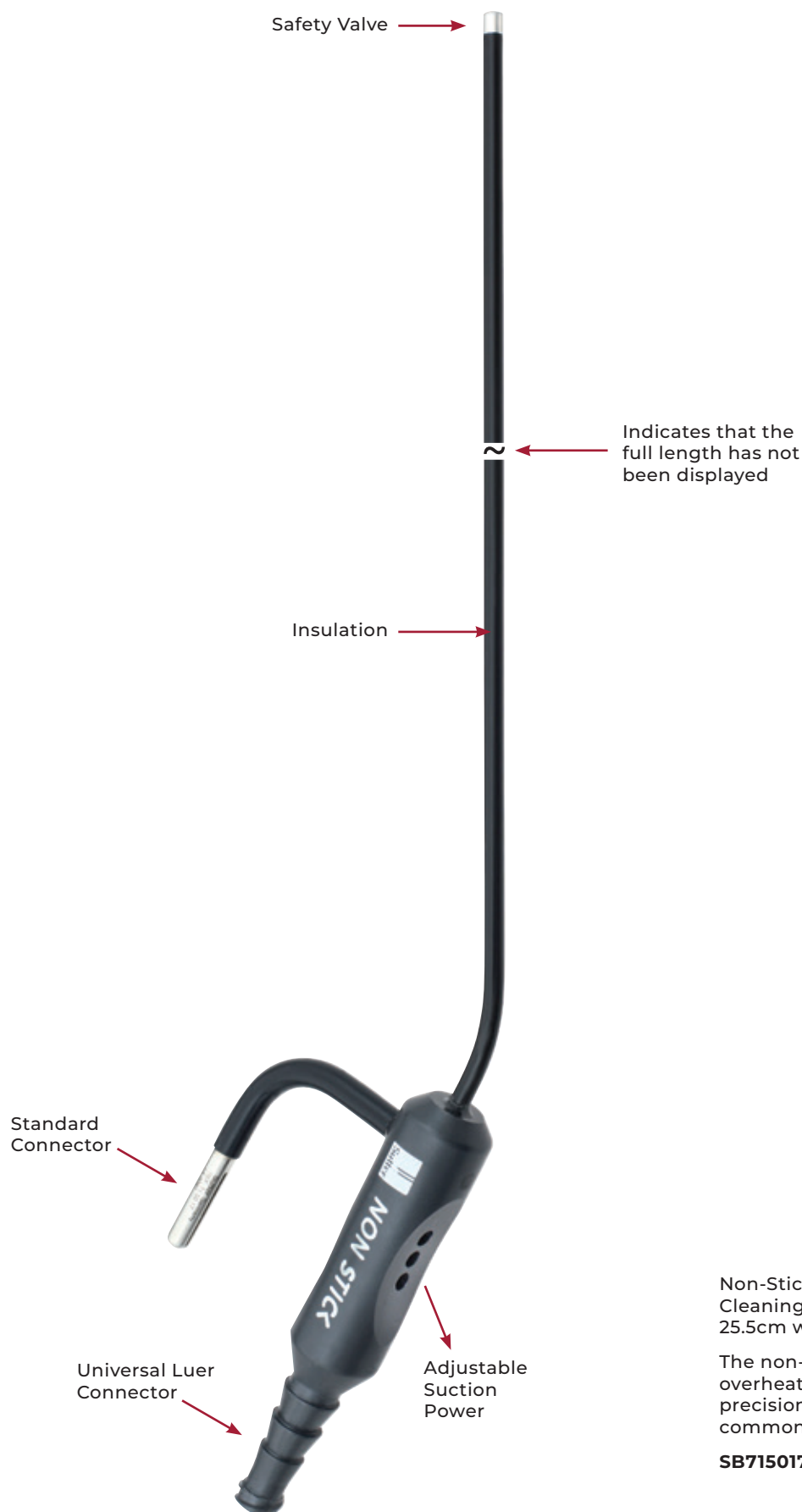


Non-Stick Monopolar Suction Tube
Cleaning brush included, $\varnothing$ 3.3mm, $\varnothing$ 2mm lumen,
25.5cm working length

The non-stick technology prevents the tips from overheating during coagulation. It is autoclavable, precision without interruptions and fits all commonly known electrosurgical units.

SB715017


SUTTER calvian® bipolar forceps

Instruments are to scale

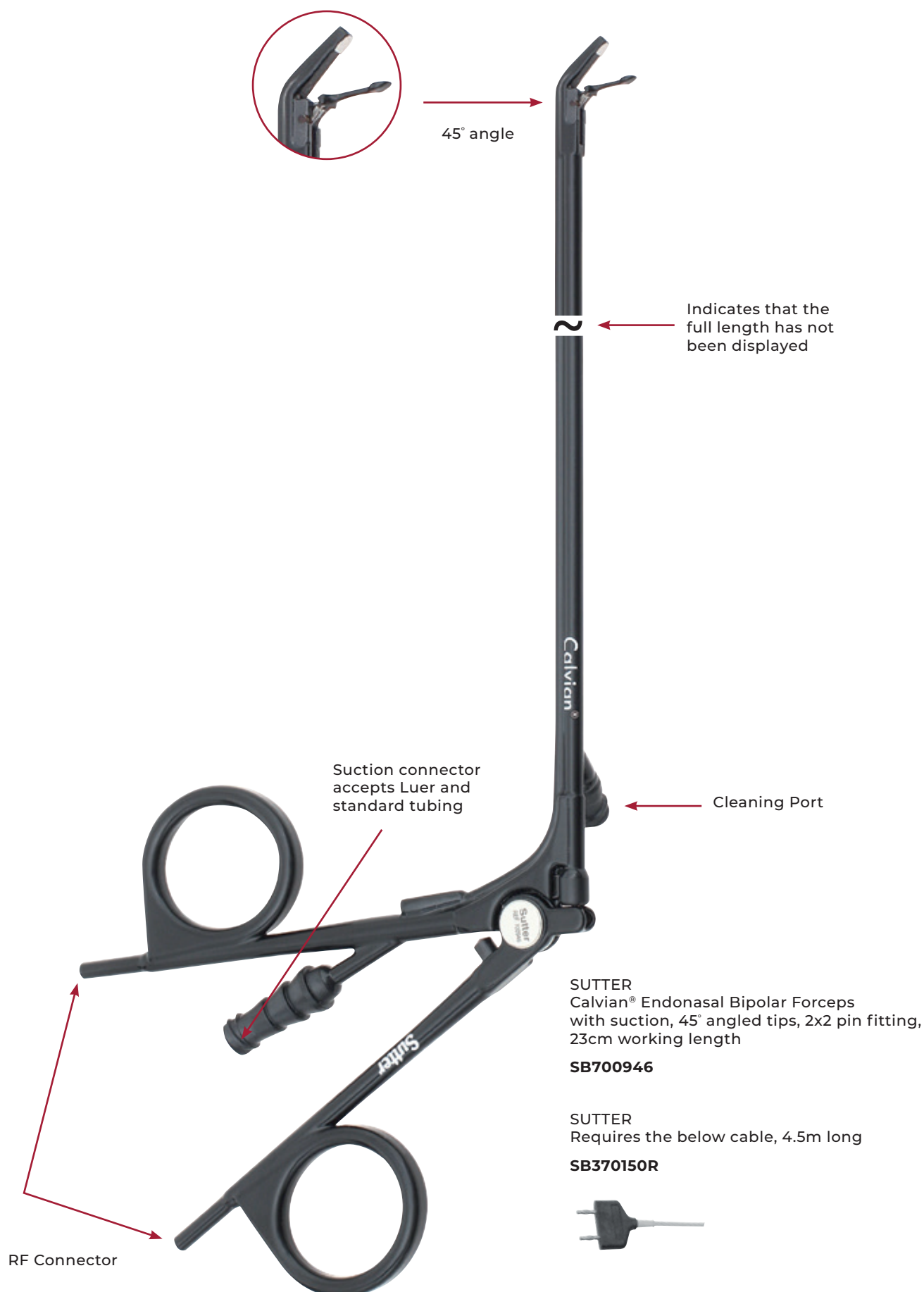

SUTTER calvian® bipolar forceps

Instruments are to scale

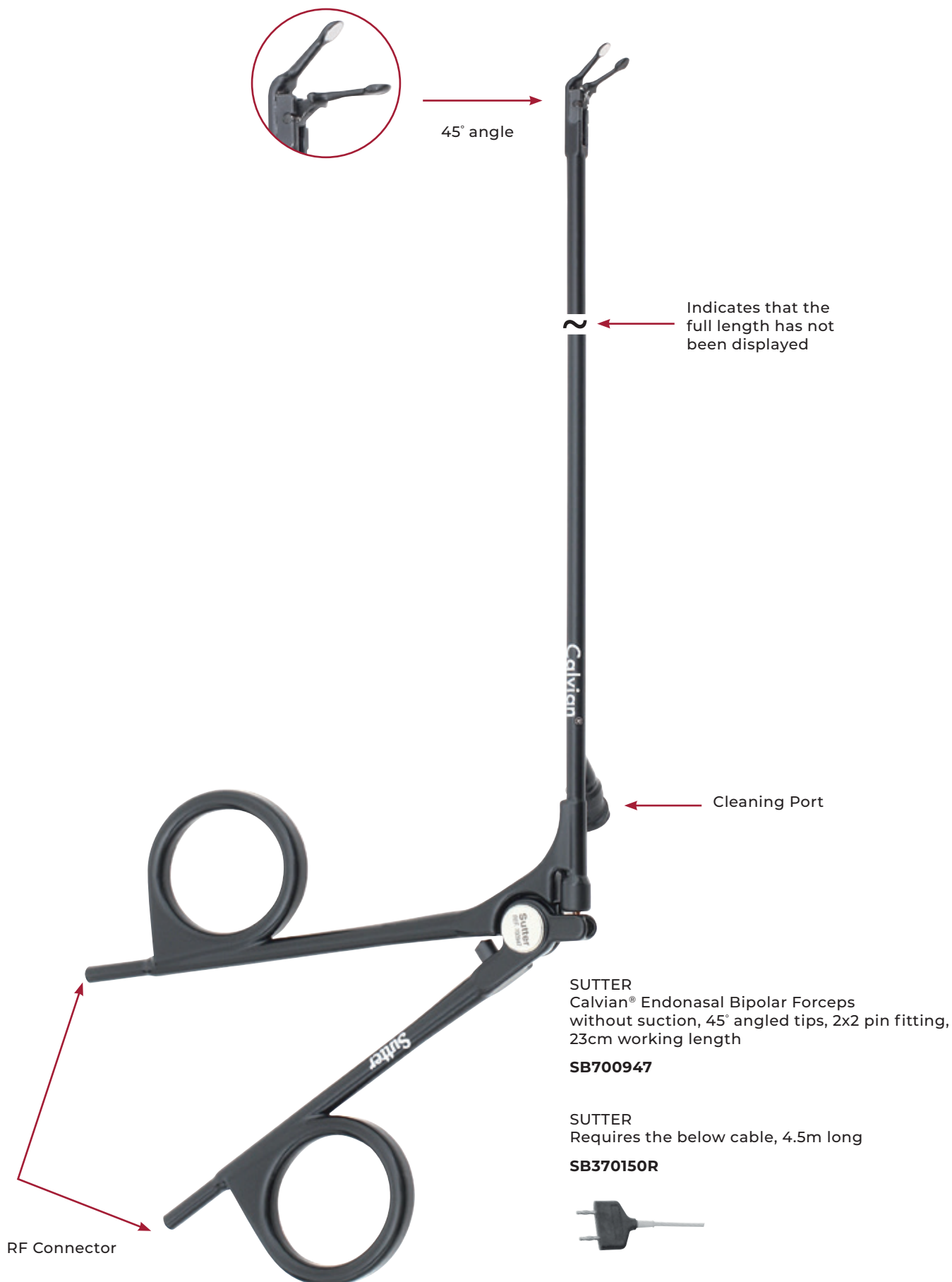

SUTTER curis® rf system

Instruments are not to scale



CURIS® - One Device, Many Applications

Meet CURIS®, our 4MHz Radiofrequency Generator for microsurgical and minimally invasive electrosurgery.



SB870010 - SUTTER Curis® Basic Set

Qty	Code	Description
1	SB360100-01	SUTTER CURIS® 4 MHz radiofrequency generator (includes mains cord, user's manual and test protocol)
1	SB360110	SUTTER footswitch two pedals for CURIS®, 4m cable
1	SB370154L	SUTTER bipolar cable for Curis®, 3m
1	SB360704	SUTTER monopolar handpiece (pencil) cut and coag, shaft Ø 2.4mm, 3m cable
1	SB360238	SUTTER cable for single-use patient plates, 3m
1	SB360222H	SUTTER safety patient plates, single-use, 50pk

SUTTER curis[®] rf system

Instruments are not to scale



SUTTER Fuego Trolley has a solid design and guarantees that the CURIS[®] 4 MHz radiofrequency generator will not shift.

Trolley comes with a hook to mount the footswitch and two storage baskets for accessories and documentation.

SB360900

SUTTER curis® for larynx surgery

Instruments are not to scale

210mm



ARROWtip™ monopolar
microdissection electrode
Basterra, autoclavable

 **A SB360371** Straight

 **D SB360374** 90°

 **B SB360372** 45°

 **E SB360375** 45°

 **C SB360373** 90°

Radiofrequency set for Larynx Surgery


	Qty	Code	Description
A	1	SB360371	ARROWtip™ monopolar microdissection electrode for larynx, straight, 236mm, 2pk
B	1	SB360372	ARROWtip™ monopolar microdissection electrode for larynx, 45° down, 233mm, 2pk
C	1	SB360373	ARROWtip™ monopolar microdissection electrode for larynx, 90° down, 231mm, 2pk
D	1	SB360374	ARROWtip™ monopolar microdissection electrode for larynx, 90° up, 234mm, 2pk
E	1	SB360375	ARROWtip™ monopolar microdissection electrode for larynx, 45° up, 236mm, 2pk

SUTTER curis[®] technology

What is Curis[®]?

The CURIS[®] 4 MHz radiofrequency generator relies on innovative 4 MHz technology. It is gentle to the tissue and effective for coagulation, for submucosal shrinkage, and for cutting. Scientific studies have shown that tissue trauma may be reduced by using CURIS[®].


Muehlifay G. et al., A study on the type of lesions archived by three electrosurgical methods and their way of heating. Romanian Journal of Morphology & Embryology, 2015, 56 (4): 1383 - 1388.

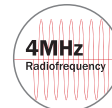

About the technology


The higher the frequency, the less the resistance of biological tissue to electromagnetic fields - up to the point where cell membranes are capacitively coupled. This effect is created by the CURIS[®] in all monopolar and bipolar modes. When using conventional electrosurgical units the magnetic field concentrates between the cells and only heats up the outer layer, however, with the CURIS[®] cell membranes are conductive, and energy is absorbed evenly inside the cells. As a result, energy is administered gently and in a highly focused fashion. Precise monopolar cuts are possible while lateral heat damage is kept to a minimum*.

*Hoffmann T.K. et al., Comparative analysis of resection tools suited for transoral robot-assisted surgery, European Archives Oto-Rhino-Laryngology. 2014, 271 (5): 1207 - 1213

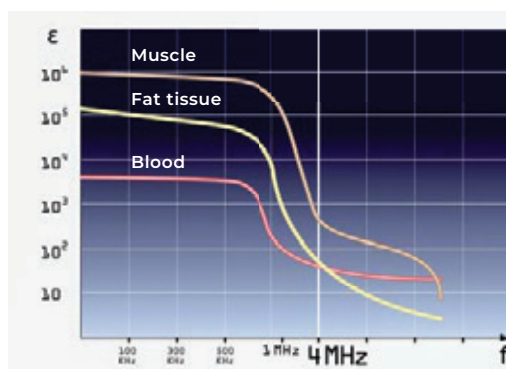


Permittivity / Frequency

This diagram shows the permittivity of tissue, which depends on the frequency of the electromagnetic field.

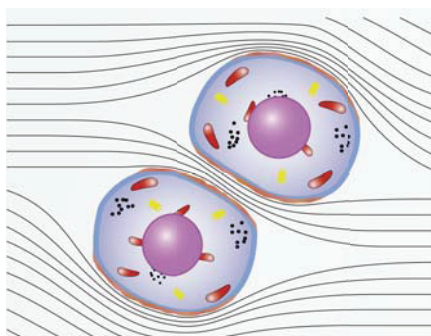

Conventional Electrosurgical Units

This electromagnetic field concentrates between the cells and heats up only the outer layer.

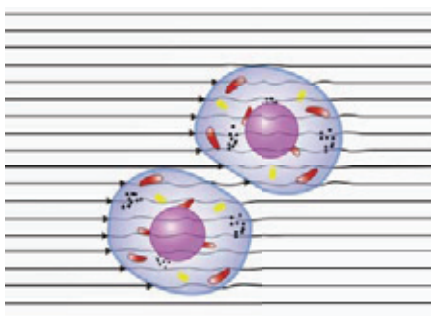

CURIS[®] 4 MHz Radiofrequency Generator

Cell membranes are conductive and the energy is absorbed evenly inside the cells. The results are highly focused tissue effects.

SUTTER curis[®] technology

AutoRF[™]


AutoRF[™] is a smart impedance control function that will tailor the power output of the CURIS[®] 4 MHz radiofrequency generator to the tissue condition. Whether it is cutting through different types of tissue (such as mucosa, muscle, fat or connective tissue) or altering tissue conditions during coagulation, the AutoRF[™] feature will deliver adapted power output as required by the different tissue impedance.

When dissecting different types of tissue in one cut (skin, fat muscles), the unit has to process and respond to the AutoRF[™] data in a flash. For this reason, the CURIS[®] has two microprocessors for additional safety and spread.



p³[™] - Technology

p³[™] stands for pulsed power performance. It is active in all coagulation modes of the CURIS[®] 4 MHz radiofrequency generator. Radiofrequency energy is delivered in about 50 small packages per second. Due to the pulsed power output, there are short breaks between the individual packages, giving the tissue enough time to absorb the energy. Highly focused, yet gentle coagulation with minimal thermal damage is possible.

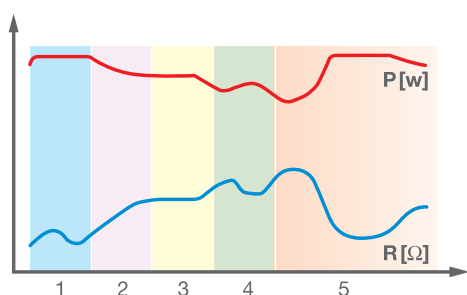

CURIS[®] Cut - Monopolar Cutting

Sections 1 to 5 show the different kinds of tissues and cutting speeds to which the unit adjusts its power output automatically.



Goldenace P/L ACN 070 452 943
ABN 28 318 606 161

Head Office
116 Cutler Road, Jandakot WA 6164

Ph 08 9417 9545
Freecall 1800 756 757

quotes@inkasurgical.com.au
www.inkasurgical.com.au